

This is an Open Access document downloaded from ORCA, Cardiff University's institutional repository:<https://orca.cardiff.ac.uk/id/eprint/99605/>

This is the author's version of a work that was submitted to / accepted for publication.

Citation for final published version:

Morgan, James , Tribble, James, Fergusson, James, White, Nicholas and Erchova, Irina 2017. The optical detection of retinal ganglion cell damage. *Eye* 31 (2) , pp. 199-205. 10.1038/eye.2016.290

Publishers page: <http://dx.doi.org/10.1038/eye.2016.290>

Please note:

Changes made as a result of publishing processes such as copy-editing, formatting and page numbers may not be reflected in this version. For the definitive version of this publication, please refer to the published source. You are advised to consult the publisher's version if you wish to cite this paper.

This version is being made available in accordance with publisher policies. See <http://orca.cf.ac.uk/policies.html> for usage policies. Copyright and moral rights for publications made available in ORCA are retained by the copyright holders.



The optical detection of retinal ganglion cell damage

JE Morgan

J Tribble

J Ferguson

N White

I Erchova

Address for Correspondence

James E Morgan
School of Optometry and Vision Sciences
Cardiff University
Maindy Road
Cardiff University
CF24 4LU

Email: morganje3@cardiff.ac.uk
Tel: 029 20874374

Word Count:

Figures 7

(Journal: Eye, Cambridge Ophthalmic Symposium)

Abstract

We provide an overview of developments in the use OCT imaging for the detection of retinal ganglion cell damage in vivo that avoid use of any exogenous ligands to label cells. The method employs high resolution OCT using broad spectral light sources to deliver axial resolution of under 5 microns. The resolution approximates that of cellular organelles, which undergo degenerative changes that progress to apoptosis as a result of axon damage. These degenerative changes are manifest as the loss of retinal ganglion cell dendrites and fragmentation of the subcellular network of organelles, in particular, the mitochondria that support dendritic structure.

These changes can alter the light scattering behaviour of degenerating neurons. Using OCT imaging techniques to identify these signals in cultured neurons, we have demonstrated changes in cultured cells and in retinal explants. Pilot studies in human glaucoma suggest that similar changes are detectable in the clinical setting.

High resolution OCT can be used to detect optical scatter signals that derive from the retinal ganglion cell layer and are associated with neuronal damage. These findings suggest that OCT instruments can be used to derive quantitative measurements of retinal ganglion cell damage. Critically, these signals can be detected at an early stage of retinal ganglion cell degeneration when cells could be protected or remodeled to support visual recovery.

Funding BBSRC UK, (BB/E017754/1), Fight for Sight (UK)

Introduction

Our understanding of retinal ganglion cell (RGC) death in diseases such as glaucoma and optic neuropathy has been transformed in recent years¹. For the most part, RGC death occurs by apoptosis², which is the culmination of a series of active cellular events involving organelle disruption. As with most chronic neuronal degenerative conditions these changes are preceded by a prolonged period of neuronal damage. In RGCs this is manifest as pruning of the processes (dendrites) that convey signals from outer retina to the RGC³⁻⁶. Dendrites process and integrate signals from bipolar cells to provide the physiological basis for the receptive field of the retinal ganglion cell- the areas in the visual field, which when stimulated will activate a single RGC. The relationship between the integrity of the receptive field and the retinal ganglion cell dendritic tree has been explored in a number of studies⁷.

Animal models of glaucoma or optic neuropathy have shown a similar pattern of RGC degeneration in which peripheral dendrites are preferentially damaged. Dendrites that are closer to the cell body (primary and secondary dendrites) are preserved in late stages of disease and^{3,8} remain discernible in eyes with end stage RGC damage{ Naskar:2003ve}. While the majority of data showing these changes are based on rodent models of disease they are also seen in other glaucoma models such as cat⁹ and primate³. Physiological studies in primates have shown that these early dendrite pruning are associated with reduced retinal ganglion cell sensitivity and contrast resolutions¹⁰.

The presence of diseased but viable cells presents an attractive target for detecting early RGC damage¹¹. RGC dendrites provide a particularly sensitive readout of neuronal degeneration in cell and tissue based studies. Recent work in explant models of RGC degeneration has highlighted the possibility that these cells can be targeted to re-establish connections with the outer retina to provide a neural substrate for the restoration of vision in diseases such as glaucoma^{11,12}. The challenge is to determine whether these changes can be discerned using non-invasive techniques that are suitable for clinical application. This paper presents data from a number of studies based on OCT imaging, which suggest that this may now be possible.

Clinical imaging for the detection of retinal ganglion cell damage

Imaging of the retinal and optic nerve has been transformed by the development of, laser based techniques that provide fast acquisition, high-resolution images of the eye. Two methods are of particular interest. Fluorescence imaging has shown particular promise for the in vivo detection of RGC degeneration^{13,14} but this remains restricted to the use of rodent models in which RGCs have been induced to fluoresce by genetic engineering. A more practical application of fluorescence imaging for the detection of RGC damage has used fluorescently tagged Annexin V antibody fragment to label cells that are in the early stages of the apoptotic pathway^{15,16}. The so-called DARC (Detection of Apoptosing Retinal Ganglion Cells) technique has achieved notable success in animal models in which RGC apoptosis has been induced by axon crush by showing a clear increase in the numbers of tagged cells following axotomy¹⁷. The technique is currently under evaluation to determine the rate of RGC apoptosis in clinical studies. The clinical utility of DARC will depend, to some extent, on precise quantification of the number of remaining RGCs. In experimental studies, the Annexin tag has been delivered by intravitreal injection¹⁷; for clinical evaluation a topical delivery method would be ideal.

Coherence based imaging

Optical coherence tomography (OCT) has become a staple clinical tool to image retinal layers in healthy and diseased tissue. The resolution of OCT instruments now falls within the scale of degenerating cellular organelles ($<5\mu\text{m}$)¹⁸ and possibly neuronal activity¹⁹. In the pre-apoptotic stages (Figure 1), mitochondrial networks fragment, preceding the nuclear disintegration which marks an irreversible commitment to cell death²⁰. While nuclear fragmentation²¹ generates significant light scatter it could only be used as an index of cell death rather than disease. By contrast, mitochondrial fragmentation is a conveniently early feature of apoptosis which could be used to highlight regions where cells are diseased but not fully committed to cell death. Work in other areas- for example scatter dependent cell sorting has suggested that changes on this scale can generate detectable levels of optical scatter²². Low coherence imaging methods have also been developed for the in vivo detection of dysplastic cellular changes most likely as a result of changes to the light scattering properties of cell nuclei²³. We have previously reported mitochondrial changes in the dendritic trees of mice with experimental (OPA1) related optic neuropathy²⁴. Early analysis of change in the mitochondrial distributions in the dendrites of midget cells in human glaucoma suggests that fragmentation and aggregation may offer feasible optical targets to generate changes in optical scatter that could be detected using OCT (Figure 2). However, the mitochondria within these dendrites will alter their distribution as a function of disease²⁵ to generate changes in optical (light) scatter that should fall within the resolution limit of high resolution OCTs. It could reasonably be argued that the availability of such high resolutions with next generation OCT devices would allow the delineation (and counting) of individual retinal ganglion cells which would be of greater clinical value. Unfortunately, the contrast that demarcates cell boundaries

requires an order of magnitude increase in sensitivity and clinical studies would be severely impacted by eye movement artifacts.

Resolution limits and definitions

For a proof of principle, the application of these techniques in the laboratory setting is simplified because exposure times can be prolonged, sample movement minimized and the aperture and optics of the test system finely controlled. A key advantage of OCT imaging is that it decouples axial and transverse resolutions²⁶. In confocal imaging systems both are hostage to the aperture and aberrations of the eye. By contrast, with OCT, the axial resolution is limited by spectral spread of the light source. Although the terms 'laser' is often used to describe the illuminating beam in an OCT, the light sources are super luminescent diodes, which deliver a spread of wavelength, approximately centred on a midpoint (mean) wavelength. Clinically used OCTs typically use light with a mean wavelength in the range 800-1040nm with a spectral spread of 50-80nm. The axial resolution is inversely related to the spectral spread with a theoretical axial resolution of 5-10 μm for these parameters.

$$l_c = \frac{2 \ln 2}{\pi} \frac{\bar{\lambda}^2}{\Delta\lambda}$$

From (²⁶), ($\bar{\lambda}$, mean wavelength, $\Delta\lambda$ the spectral width).

It can be seen that an increase in the spectral spread ($\Delta\lambda$), will reduce the minimum resolvable dimension (l_c), resulting in an increase in the resolution. For example, OCTs used for experimental work (e.g. Phoenix Research Laboratories) have a spectral spread of the order of 140nm, to provide a theoretical axial resolution of approximately 3 μm). Until recently, cost has been a limiting factor for the use of wider spectral light sources (e.g. Titanium Sapphire lasers). The wider availability of multiplexed diodes whose wavelength output can be approximated to a Gaussian distribution with a large spectral spread (> 100 nm) will help in bringing higher resolution OCT into the clinical domain. When the resolution is matched in transverse and axial directions it is called isotropic- the ideal state for an imaging device. For most OCTs the axial resolution exceeds the transverse resolution since the latter is affected by aberrations of the human eye. These can be overcome using adaptive optics to yield high-resolution isotropic resolutions²⁷. These OCT systems are however, complex and with limited views of the retina and optic nerve which can limit clinical application.

Light scatter during neuronal degeneration

Since OCTs can potentially detect signals from objects approximating the dimension of subcellular organelles, we reasoned that they could be used to detect signals that correlate with the time course of neuronal degeneration. OCT images contain noise (speckle) arising from use of a laser as a light source but also from interference with multiple back scattering objects in a region of interest. The sources of speckle in OCT imaging has been covered in depth by Schmitt et al.²⁸ who usefully differentiated speckle into that containing information

relevant to the imaged objects from speckle that degraded the image quality (Figure 3a). We therefore developed methods to extract useful measurements of optical texture from the speckle associated with OCT images. The signal to support this analysis combines optical scatter from changes in the tissue under examination, speckle due to forward light scatter and noise in the imaging system and detector. These signals are usually removed by image processing in commercial OCTs to facilitate segmentation of retinal boundaries.

We therefore constructed a custom OCT coupled to a standard compound microscope to examine changes in texture that occur in both cultured retinal ganglion cells and retinal explants²⁹ (Figure 3b). The OCT used a Femtosecond Integral Ti:Sapphire laser with a high bandwidth of 140 nm to provide a theoretical bandwidth of 3 μ m. The examination of single cells over time allowed the detection of subcellular changes that occur once cell death has been initiated. We selected RGC5 cells which at the time of the study were thought to provide a useful mimic of retinal ganglion cells. Subsequent analysis has shown that these are most likely a murine neuronal cell line which can express some markers shared by retinal ganglion cells when suitably induced³⁰.

RGC5 neuronal degeneration over time

Histological examination confirmed mitochondrial aggregation and caspase activation in RGC5 cell death following the initiation of apoptosis by the addition of staurosporine^{29,31}. The details of the mathematical modeling are described in details elsewhere²⁹. In brief, the image was parameterized using 65 features, which were then selected on the basis of the accuracy with which they correlated with culture stage. The selection of these dimensions was hypothesis free and based on the (unprocessed) spectral signal from the OCT. The aim of these pilot experiments was to determine if texture features could be detected that correlated with the stage of degeneration in the RGCs and allowed the effective discrimination of healthy from dying cells. Images were taken at regular intervals of the same region of interest, prior to RGC5 apoptosis (as demonstrated by Tunel staining).

We employed a Gaussian Mixture Model analysis to identify those variables related to image texture which corresponded the degree of cellular degeneration. We aggregated these measures into a single dimensionless coefficient – the Texture Index (TI) to plot against time following the administration of staurosporine. These parameters were then fixed and used to analyse separate populations of RGC5 cells following the administration of staurosporine and agent that can induce apoptosis. We observed a progressive shift in the TI within 20 minutes of staurosporine, prior to the onset of significant apoptosis (Figure 4). These observations were consistent with a change in light scatter as a function of neuronal degeneration.

We next examined whether similar signals could be detected in in retinal explants. Explants provide an excellent model for the study of retinal ganglion cell degeneration since RGC

degeneration is initiated by axotomy during specimen preparation³². We have analysed the changes in RGC dendritic structure using particle-mediated fluorescent labeling (Diolistics), showing that dendrite pruning occurs within 6 hours of axotomy. These changes are shown in figure 5 demonstrating that RGC dendrites are extensively pruned during the first 3 days of culture,¹² supporting their use as a model of dendrite degeneration. For the purposes of this study, we randomly selected regions of interest from the inner plexiform layer (IPL) and subjected these to texture analysis using similar parameter settings as used for the RGC5 study. Retinal explants were scanned enfaces, retinal ganglion cell side up, and maintained their optical clarity throughout the experiment. The TI showed a monotonic increase as a function of time consistent with observations in single cell data (Figure 6).

Human studies

Subsequent studies determined whether texture analysis could confer any benefit in the detection of glaucoma at an early stage in the disease³³. The experimental OCT device, which had a lower spectral spread and longer wavelength compared with the RGC5 study (NP photonics ASE light source, 1024nm, FWHM 80nm, theoretical axial resolution 7µm), was connected to a clinical scanner and macular scans taken of patients with early/moderate glaucoma (14 normal, 10 glaucoma). As for the explant studies, texture signatures were extracted from the inner plexiform layer but in addition these were passed to a Supported Vector Machine (SVM) methods for classification to determine if they could confer any additional diagnostic precision in the discrimination of patients with early/moderate glaucoma. Several images processing steps were required to eliminate artifacts arising from micro saccades³⁴. A small but potentially useful improvement in diagnostic accuracy (up to 4%) could be seen when texture features were added to more conventional parameters used to discriminate glaucomatous eyes.³³

Texture analysis for the detection of retinal ganglion cell disease

Our studies and the work of others support the use of texture analysis to quantify some dimensions of neuronal health taking RGCs as a specific example. The cell-based studies were performed with a high resolution OCT employing broad spectral light sources. Since the costs of these light sources have reduced substantially their use in the clinical setting is to be anticipated. All studies were conducted on custom OCTs in which all stages of image processing were known; this is not necessarily the case with commercially available devices, which could limit further clinical application of this approach.

Our data suggest that texture measurements and indices could be derived if OCT image filtering and segmentation were undertaken as suggested in figure 7. In this model, unfiltered data is stored and processed to derive texture data in one analytic pathway and filtered data

passed to another to segment retinal layers and define regions of interest. The texture scores can then be aligned to the segmented images to provide an index of neuronal health for selected retinal layers. In the context of studying RGC health we suggest the IPL as an informative region of interest. It is clear that this technique has potential for the quantification of cellular health in other retinal layers, for example the photoreceptor layer for patient at risk of macular degeneration.

Figure legends

Figure 1

Overview of the sequence of cellular changes that precede apoptotic (programmed) cell death. 'Disease' covers changes that are reversible. 'Death' covers the point at which the cells commit to cell death. DARC: detection of apoptosing retinal ganglion cells.

Figure 2

Serial block face scanning electron microscope reconstruction images of human midget RGC from a normal and glaucomatous eye. The processes in the glaucomatous eye are reduced and mitochondria show aggregation within the remaining dendrite structure. Each RGC serially sectioned at 100nm. Scale 10 μ m.

Figure 3a

Schematic of the OCT setup. Light source mean wavelength 800nm (bandwidth 140 nm). From ²⁹

Figure 3b

Origin of scatter (speckle signal) shown the interaction of forward and back scattering of light to generate speckle. Adapted from ²⁸

Figure 4

Gaussian Mixture Model analysis of early stages of apoptosis. Cell classification was based at each time point on the basis of texture analyses. A progressive shift in the Texture Index (TI) defined as before was observed within 20 min of the administration of staurosporine. Significance levels were calculated by t-test (*p,0.05, ** p,0.01,*** p,0.001). From ²⁹

Figure 5

Sholl plots representing the complexity of RGC dendritic trees as a function of time in culture for murine retinal explants. The arrow indicates the leftward shift in the distributions as a function of time. Note that the first changes are seen within 6 hours of culture. From ¹²

Figure 6

Changes in texture index as a function of time for murine retinas maintained in culture. OCT measurements were taken from the IPL within 1 hour of explant preparation. The texture index shows a steady increase with time, changing as early as 10 minutes of culture time. The inset diagram shows the orientation of the OCT scanner and retinal explant. From ²⁹

Figure 7

Scheme showing how OCT data can be separated into filtered data for the detection of retinal boundaries and unfiltered data for texture analysis. The two analyses can be merged to generate layer specific texture maps.

References

1. Almasieh M, Wilson AM, Morquette B, Cueva Vargas JL, Di Polo A. The molecular basis of retinal ganglion cell death in glaucoma. *Progress in retinal and eye research*. 2012;31(2):152-181.
2. Kerrigan LA, Zack DJ, Quigley HA, Smith SD, Pease ME. TUNEL-positive ganglion cells in human primary open-angle glaucoma. *Archives of Ophthalmology*. 1997;115(8):1031-1035.
3. Weber AJ, Kaufman PL, Hubbard WC. Morphology of single ganglion cells in the glaucomatous primate retina. *Invest Ophthalmol Vis Sci*. 1998;39(12):2304-2320.
4. Morgan JE, Datta AV, Erichsen JT, Albon J, Boulton ME. Retinal ganglion cell remodelling in experimental glaucoma. *Adv Exp Med Biol*. 2006; 572:397-402.
5. Kalesnykas G, Oglesby EN, Zack DJ, et al. Retinal ganglion cell morphology after optic nerve crush and experimental glaucoma. *Invest Ophthalmol Vis Sci*. 2012;53(7):3847-3857.
6. Williams PA, Howell GR, Barbay JM, et al. Retinal ganglion cell dendritic atrophy in DBA/2J glaucoma. *PloS one*. 2013;8(8):e72282. doi:10.1371/journal.pone.0072282.
7. Brown SP, He S, Masland RH. Receptive field microstructure and dendritic geometry of retinal ganglion cells. *Neuron*. 2000;27(2):371-383.
8. Williams PA, Morgan JE, Votruba M. Opa1 deficiency in a mouse model of dominant optic atrophy leads to retinal ganglion cell dendropathy. *Brain*. 2010;133(10):2942-2951.
9. Shou T, Liu J, Wang W, ZHOU Y. Differential dendritic shrinkage of α and β retinal ganglion cells in cats with chronic glaucoma. *Invest Ophthalmol Vis Sci* 2003;44(7) 3005-3010.
10. Weber AJ, Harman CD. Structure-function relations of parasol cells in the normal and glaucomatous primate retina. *Invest Ophthalmol Vis Sci*. 2005;46(9):3197-3207.
11. Morgan JE. Retina ganglion cell degeneration in glaucoma: an opportunity missed? A review. *Clin Experiment Ophthalmol*. 2012;40(4):364-368.
12. Binley KE, Ng WS, Barde Y-A, Song B, Morgan JE. Brain-derived neurotrophic factor prevents dendritic retraction of adult mouse retinal ganglion cells. Kirik D, ed. *European Journal of Neuroscience*. 2016;44(3):2028-2039. doi:10.1111/ejn.13295.
13. Leung CK-S, Weinreb RN, Li ZW, et al. Long-term in vivo imaging and measurement of dendritic shrinkage of retinal ganglion cells. *Invest Ophthalmol Vis Sci*. 2011;52(3):1539-1547.
14. Leung CK-S, Lindsey JD, Crowston JG, Lijia C, Chiang S, Weinreb RN. Longitudinal profile of retinal ganglion cell damage after optic nerve crush with blue-light confocal scanning laser ophthalmoscopy. *Invest Ophthalmol Vis Sci*. 2008;49(11):4898-4902.
15. Cordeiro MF, Guo L, Luong V, et al. Real-time imaging of single nerve cell apoptosis in retinal neurodegeneration. *Proc Natl Acad Sci USA*. 2004;101(36):13352-13356.
16. Cordeiro M. DARC: a new method for detecting progressive neuronal death. *Eye*.

- 2007 21 S15-17.
17. Guo L, Cordeiro MF. Assessment of neuroprotection in the retina with DARC. In: *Glaucoma: an Open Window to Neurodegeneration and Neuroprotection*. Vol 173. Progress in Brain Research. Elsevier; 2008:437-450.
 18. Drexler W. State-of-the-art retinal optical coherence tomography. *Progress in retinal and eye research*. 2008 27: 45-88.
 19. Bizheva K, Pflug R, Hermann B, et al. Optophysiology: depth-resolved probing of retinal physiology with functional ultrahigh-resolution optical coherence tomography. *Proc Natl Acad Sci USA*. 2006;103(13):5066-5071.
 20. Arnoult D. Mitochondrial fragmentation in apoptosis. *Trends in Cell Biology*. 2007;17(1):6-12.
 21. Farhat G, Mariampillai A, Yang VXD, Czarnota GJ, Kolios MC. Optical coherence tomography speckle decorrelation for detecting cell death. Wax AP, Backman V, eds. SPIE 2011;7907: doi:10.1117/12.874682.
 22. Darzynkiewicz Z, Bruno S, Del Bino G, et al. Features of apoptotic cells measured by flow cytometry. *Cytometry Part A*. 1992;13(8):795-808.
 23. Wax A, Chalut KJ. Nuclear morphology measurements with angle-resolved low coherence interferometry for application to cell biology and early cancer detection. *Anal Cell Pathol* 2011. 34(5) 207-222..
 24. Williams PA, Piechota M, Ruhland von C, Taylor E, Morgan JE, Votruba M. Opa1 is essential for retinal ganglion cell synaptic architecture and connectivity. *Brain*. 2012;135(2):493-505.
 25. Mizielińska S, Frenguelli B, Harvey J. Mitochondrial dysfunction and dendritic beading during neuronal toxicity. *J. Biol. Chem*. 282(36) 26235-26244.
 26. Fercher AF, Drexler W, Hitzenberger CK, Lasser T. Optical coherence tomography - principles and applications. *Rep Prog Phys*. 2003;66(2):239-303.
 27. Torti C, Povazay B, Hofer B, et al. Adaptive optics optical coherence tomography at 120,000 depth scans/s for non-invasive cellular phenotyping of the living human retina. *Opt Express*. 2009;17(22):19382-19400.
 28. Schmitt JM, Xiang SH, Yung KM. Speckle in Optical Coherence Tomography. *J Biomed Opt*. 1999;4(1):95-105. doi:10.1117/1.429925.
 29. Tudor D, Kajić V, Rey S, et al. Non-Invasive Detection of Early Retinal Neuronal Degeneration by Ultrahigh Resolution Optical Coherence Tomography. *PloS one*. 2014;9(4):e93916.
 30. Van Bergen NJ, Wood JPM, Chidlow G, et al. Recharacterization of the RGC-5 Retinal Ganglion Cell Line. *Invest Ophthalmol Vis Sci*. 2009;50(9):4267-4272.
 31. Kim SC, Kang TJ. Texture classification and segmentation using wavelet packet frame and Gaussian mixture model. *Pattern Recognition*. 2007;40(4):1207-1221.
 32. Johnson TV, Martin KR. Development and Characterization of an Adult Retinal Explant Organotypic Tissue Culture System as an In Vitro Intraocular Stem Cell Transplantation Model. *Invest Ophthalmol Vis Sci*. 2008;49(8):3503-3512.
 33. Anantrasirichai N, Achim A, Morgan JE, Erchova I, Nicholson L. SVM-based texture

- classification in Optical Coherence Tomography. In: IEEE; 2013:1332-1335.
34. Anantrasirichai N, Nicholson L, Morgan JE, et al. Adaptive-weighted bilateral filtering and other pre-processing techniques for optical coherence tomography. *Computerized Medical Imaging and Graphics*. 2014;38(6):526-539.

Figure 1

Cell changes

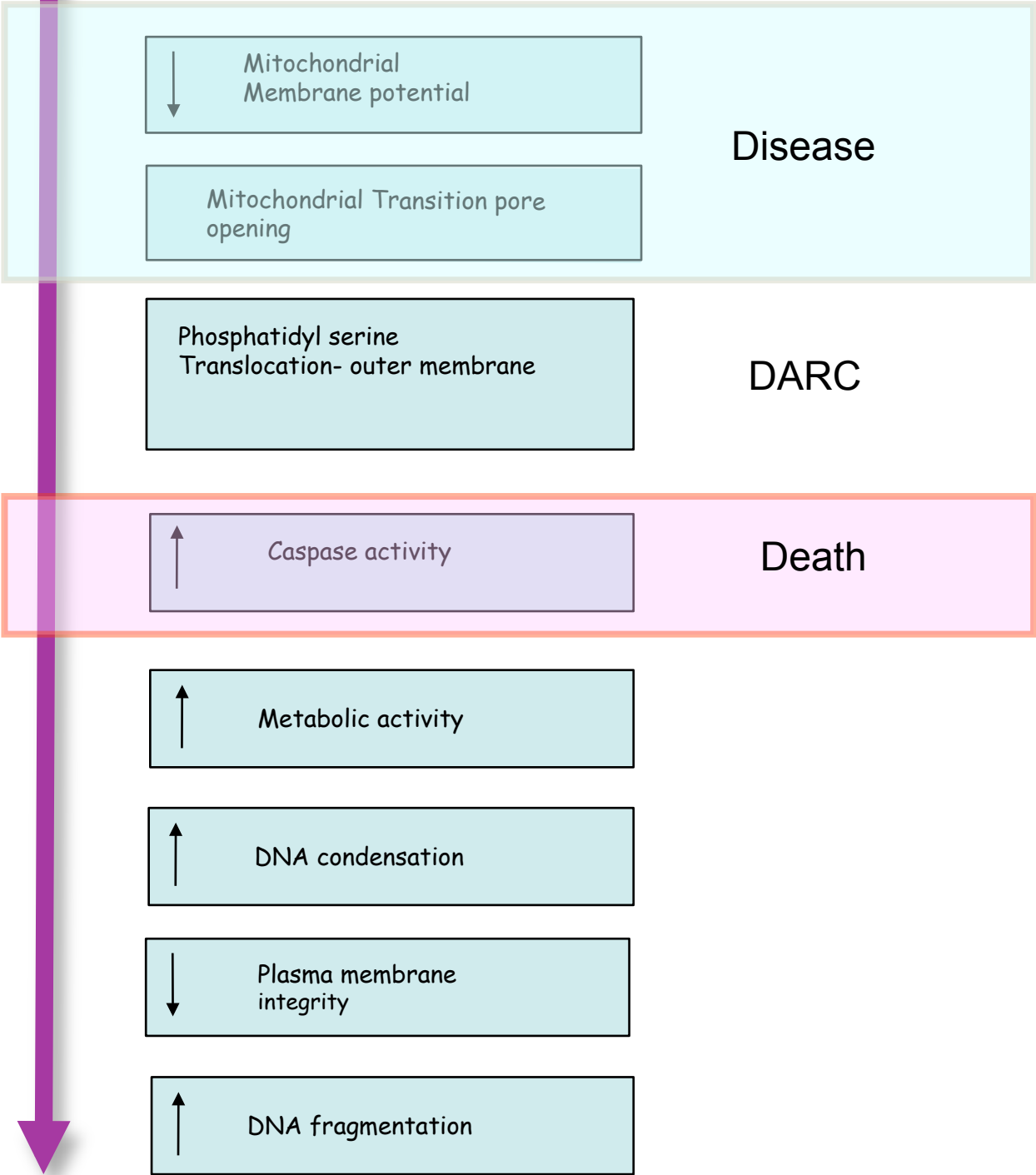


Figure 2

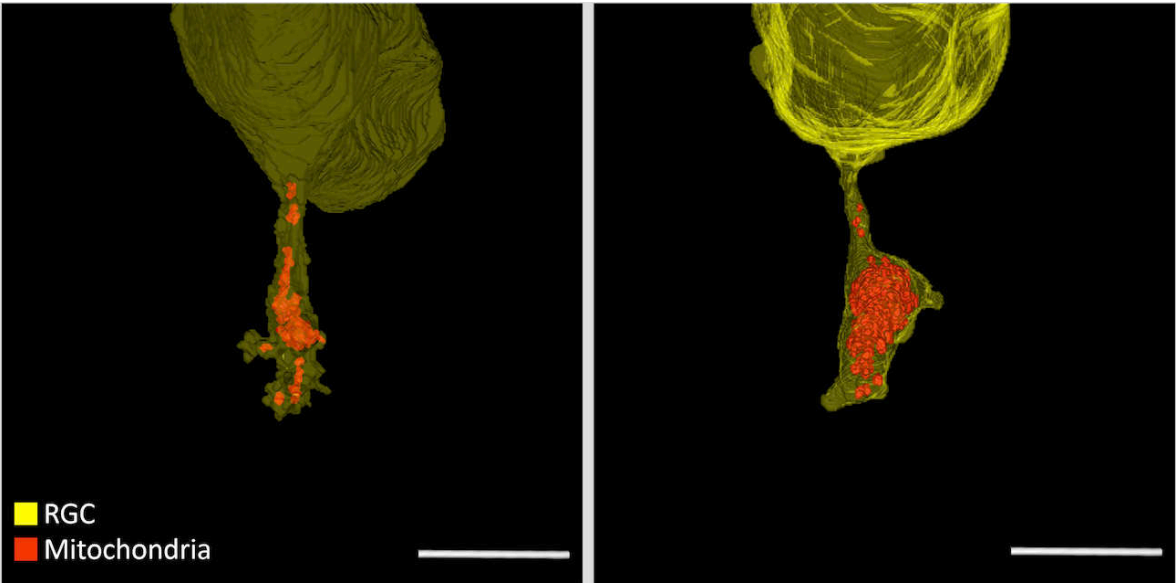


Figure 3a

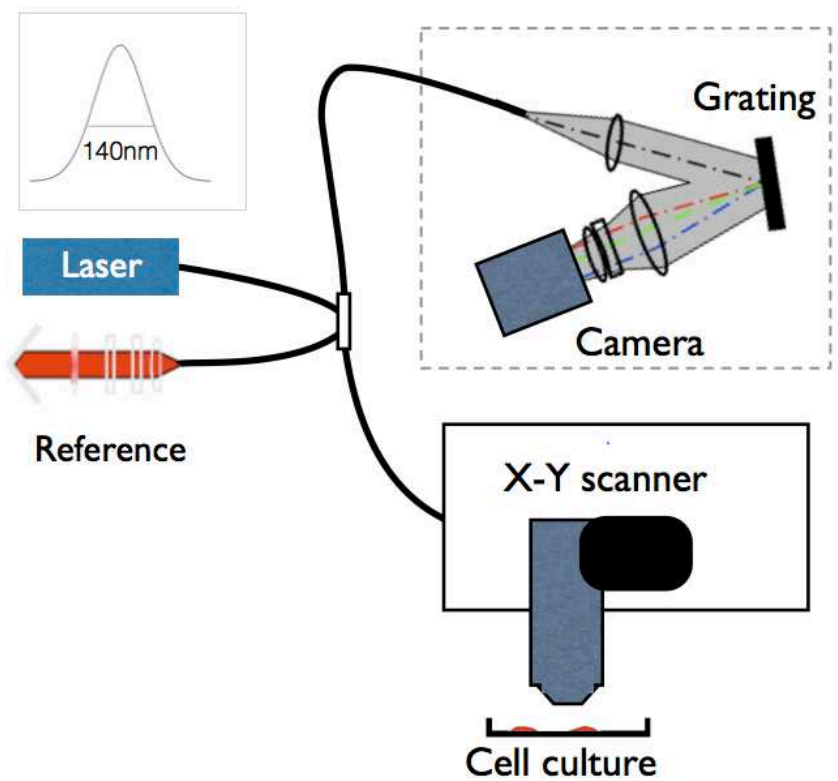


Figure 3b

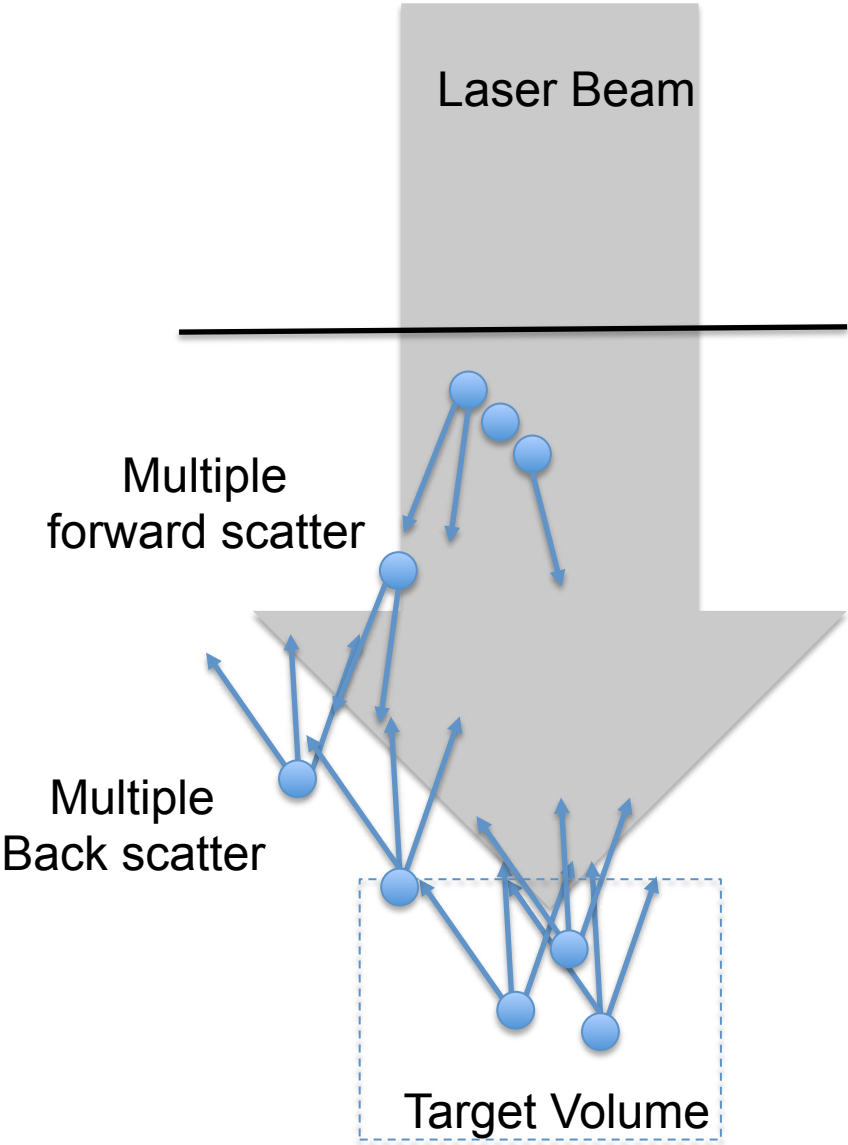


Figure 4

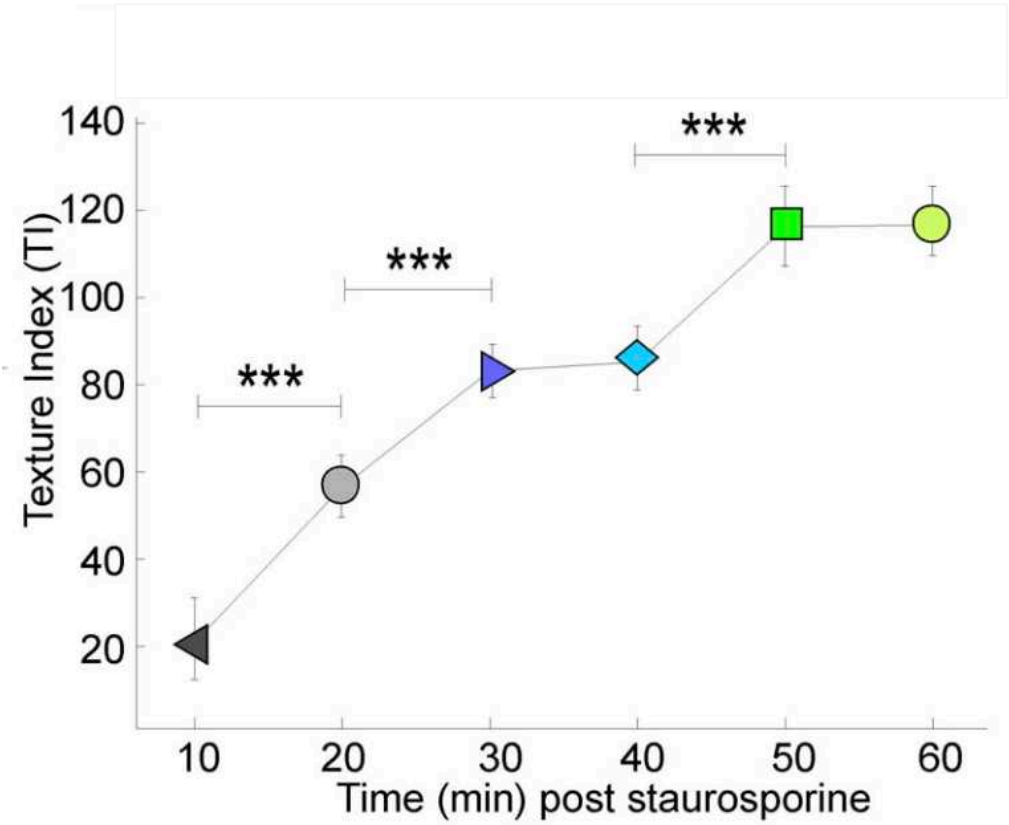


Figure 5

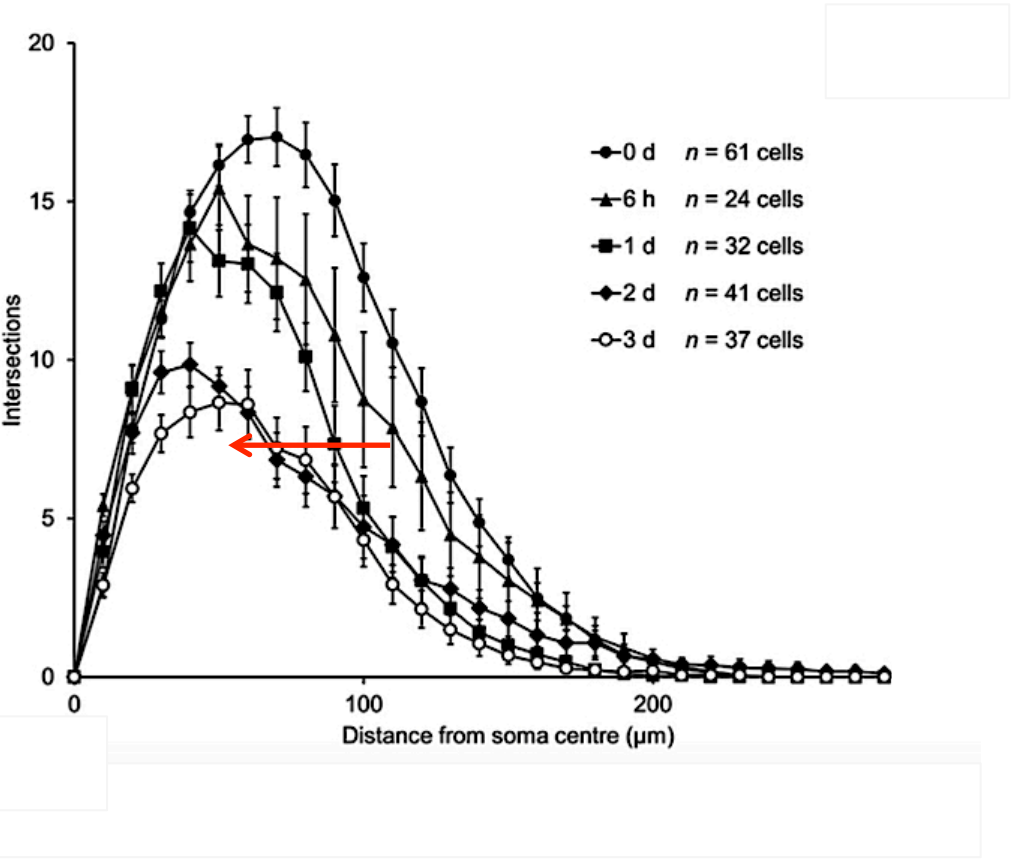


Figure 6

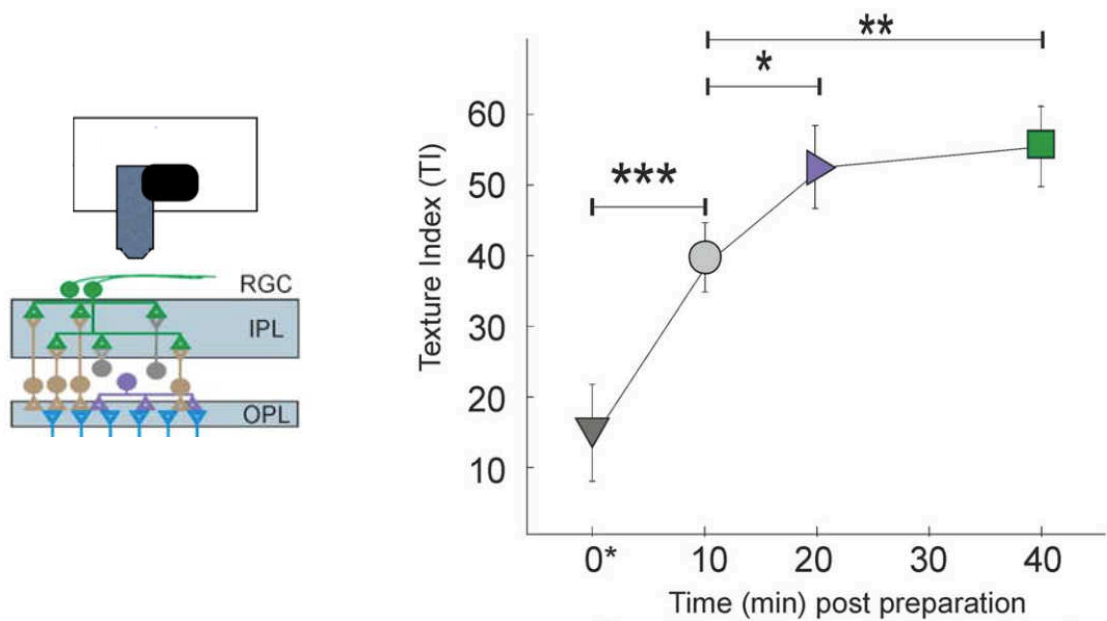


Figure 7

