

SMALL BOWEL

Blood-Based T-Cell Diagnosis of Celiac Disease



Olivia G. Moscatelli,^{1,2} Amy K. Russell,¹ Lee M. Henneken,³ Linda Fothergill,¹ Allan Motyer,⁴ Hugh Reid,⁵ Jamie Rossjohn,^{5,6} Vanessa Bryant,^{1,2} Robert P. Anderson,⁷ Melinda Y. Hardy,^{1,2} and Jason A. Tye-Din^{1,2,3}

¹Immunology Division, The Walter and Eliza Hall Institute, Parkville, Victoria, Australia; ²Department of Medical Biology, University of Melbourne, Parkville, Victoria, Australia; ³Department of Gastroenterology, The Royal Melbourne Hospital, Parkville, Victoria, Australia; ⁴Advanced Technology and Biology, The Walter and Eliza Hall Institute, Parkville, Victoria, Australia; ⁵Infection and Immunity Program and Department of Biochemistry and Molecular Biology, Biomedicine Discovery Institute, Monash University, Clayton, Victoria, Australia; ⁶Institute of Infection and Immunity, Cardiff University School of Medicine, Heath Park, Cardiff, United Kingdom; and ⁷Gastroenterology Department, Mackay Base Hospital, Mackay, Queensland, Australia

BACKGROUND & AIMS: Current diagnosis of celiac disease (CeD) is inaccurate in patients consuming a gluten-free diet (GFD). Blood-based diagnostics targeting gluten-specific T cells, such as tetramer assays, are highly sensitive and specific but are impractical for clinical use. We evaluated the potential of a simple, whole-blood assay measuring interleukin-2 (IL-2) release (WBAIL-2) for detecting gluten-specific T cells to aid in CeD diagnosis. **METHODS:** WBAIL-2 was assessed in 181 adults; 88 with CeD (75 consuming a GFD, 13 consuming gluten) and 93 controls (32 consuming a GFD with nonceliac gluten sensitivity, 61 healthy). In vitro IL-2 release in whole blood after gluten peptide stimulation was measured. The assay's performance was compared with tetramer-based methods, and serum IL-2 levels were monitored before and after a single-dose gluten challenge. Correlations between IL-2 levels, tetramer-positive T-cell frequencies, and symptoms were examined. **RESULTS:** The WBAIL-2 assay demonstrates high accuracy for CeD diagnosis, even in patients consuming a strict GFD. Optimized dual cutoffs in HLA-DQ2.5⁺ patients showed high sensitivity (90%) and specificity (95%), with lower sensitivity (56%) in HLA-DQ8⁺ CeD. WBAIL-2 correlated strongly with the frequency of tetramer-positive gluten-specific CD4⁺ T cells and serum IL-2 levels after a gluten challenge. Elevated WBAIL-2 levels predicted gluten-induced symptom severity, such as vomiting. The assay required only small blood volumes and performed comparably with tetramer-based methods. **CONCLUSIONS:** Gluten-stimulated IL-2 secretion indicates the presence of pathogenic gluten-specific CD4⁺ T cells and is a useful diagnostic for CeD. WBAIL-2 and serum IL-2 after gluten could be complementary and allow biopsy-free CeD diagnosis. WBAIL-2 may help diagnose and monitor other CD4⁺ T cell-driven diseases.

but global prevalence is high (1%–2%), with 50% to 80% of cases undiagnosed or diagnosed late.^{2–4}

Reliable diagnosis of CeD requires patients to consume gluten, because first-line serology tests and confirmatory endoscopic duodenal biopsy specimens depend on the detection of gluten-induced injury, which normalizes on a GFD. This requirement poses challenges for patients already avoiding gluten, and with the anticipated approval of nondietary therapies, accurate tools for confirming diagnosis and monitoring disease severity will become increasingly important.⁵

The diagnosis and therapeutic monitoring of CeD and other chronic immune diseases, as well as infectious diseases and cancers, could be significantly advanced by measuring disease-relevant antigen-specific effector CD4⁺ T cells. However, this approach has seen limited clinical application outside of academic studies, vaccine development, and pharmaceutical research. Tuberculosis is one of the few conditions diagnosed by measuring antigen-specific T cells. Circulating mycobacteria-specific T cells are relatively abundant in latent tuberculosis and can be detected by enzyme-linked immunospot (ELISpot) or whole-blood cytokine release assay.^{6,7}

Unfortunately, the frequency of antigen-specific T cells in blood without antigen rechallenge is below the detection limits of these platforms for most diseases, and in many cases, the identity of the triggering antigenic peptides is poorly defined. Although flow cytometry with peptide-major histocompatibility complex multimers has the sensitivity to detect these rare cells, it requires knowledge of the relevant T-cell epitopes and is restricted to highly specialized laboratories.

Keywords: Celiac Disease; T cells; Interleukin-2; Diagnostics.

Celiac disease (CeD) is caused by an adaptive immune response to dietary gluten mediated by CD4⁺ gluten-specific T cells that results in chronic small intestinal injury and diverse chronic morbidities.¹ The only available treatment is a strict, life-long, gluten-free diet (GFD). Expedient diagnosis is important to minimize long-term complications,

Abbreviations used in this paper: APC, allophycocyanin; AUC, area under the curve; CeD, celiac disease; CI, confidence interval; CLIP, class II-associated invariant chain peptide; ELISpot, enzyme-linked immunospot; CV, coefficient of variation; GCIL-2, gluten challenge interleukin-2; GFD, gluten-free diet; HLA, human leukocyte antigen; IFN, interferon; IL, interleukin; NCGS, nonceliac gluten sensitivity; PBMC, peripheral blood mononuclear cells; PE, phycoerythrin; T_{EM}, T-effector memory; TNF, tumor necrosis factor; WBA, whole-blood assay; WBAIL-2, whole-blood assay measuring interleukin-2 release.

Most current article

© 2025 The Author(s). Published by Elsevier Inc. on behalf of the AGA Institute. This is an open access article under the CC BY license (<http://creativecommons.org/licenses/by/4.0/>).

0016-5085

<https://doi.org/10.1053/j.gastro.2025.05.022>

WHAT YOU NEED TO KNOW**BACKGROUND AND CONTEXT**

Celiac disease (CeD) diagnosis is challenging, particularly for patients on gluten-free diets. A whole blood interleukin-2 (IL-2) release assay (WBAIL-2) may simplify detection of gluten-specific T cells and be a straightforward diagnostic tool for CeD.

NEW FINDINGS

WBAIL-2 is a simple and accurate CeD diagnostic even in patients consuming a gluten-free diet. IL-2 production reflects the activation of gluten-specific T cells, correlates with their frequency, and stratifies symptom severity to gluten.

LIMITATIONS

This was a single-center study with relatively small subgroup sizes and untested reproducibility across laboratories. Pediatric populations and patients on immunosuppressants were not assessed, and prospective validation is required.

CLINICAL RESEARCH RELEVANCE

WBAIL-2 is an accurate diagnostic tool for CeD, even in patients consuming a gluten-free diet. It matches the sensitivity of major histocompatibility complex-tetramer detection but requires significantly less blood and is easier to perform. It has the potential to replace gastroscopy and small intestinal biopsies. WBAIL-2 could aid in diagnosing and monitoring other diseases driven by CD4⁺ T cells.

BASIC RESEARCH RELEVANCE

This study confirms IL-2 as a sensitive and specific biomarker of gluten-specific T-cell activation, capable of detecting a single gluten-specific T cell in 4 mL of blood. WBAIL-2's practical advantages, including minimal blood requirements and no need for specialized infrastructure, make it a valuable tool for studying immune responses in CeD and other T cell-driven disorders.

CeD is an appealing candidate for a T-cell diagnostic due to the central role of gluten-specific CD4⁺ T cells in its pathogenesis and the clear identification of immunodominant gluten epitopes driving these cells.⁸ Gluten-specific CD4⁺ T cells are specific to CeD; however, they circulate in low numbers in the blood and can only be detected by peptide-major histocompatibility complex tetramers^{9,10} or by functional T-cell assays, such as ELISpot, which require a preceding oral gluten challenge.^{11,12}

Recent advances suggest that ultrasensitive assays detecting interleukin-2 (IL-2) release after gluten ingestion could complement serology and enable non-biopsy specimen-diagnosis of CeD.^{13–17} IL-2 release has demonstrated high sensitivity and specificity in distinguishing CeD patients from healthy controls and nonceliac gluten sensitivity (NCGS).^{18,19} A whole-blood assay for IL-2 (WBAIL-2) was developed, in which gluten peptides stimulate blood in vitro, and IL-2 release is measured.²⁰

This study assessed the feasibility and performance of these novel, but technically straightforward, in vivo and

in vitro gluten-stimulated IL-2 release assays in CeD. The findings provide proof of concept for these ultrasensitive assays being comparable to human leukocyte antigen (HLA)-DQ:gluten peptide tetramer assays in detecting rare antigen-specific CD4⁺ T cells and having clinical utility in the diagnosis of CeD.

Materials and Methods**Study Design and Participants**

This was a single-center, investigator-led study performed at the Walter and Eliza Hall Institute, with recruitment via the Royal Melbourne Hospital. Adult participants aged 18 to 75 years with medically diagnosed CeD (treated and active cohorts), NCGS, and healthy controls were recruited. All participants provided written informed consent. The study was approved by the Human Research Ethics Committees of the Royal Melbourne Hospital (2021.210) and the Walter and Eliza Hall Institute (21/18) and was conducted in accordance with the ethical principles in the Declaration of Helsinki.

Inclusion criteria for CeD participants was documentation confirming a past CeD diagnosis based on duodenal villous atrophy (Marsh 3), positive celiac serology, and supportive clinical criteria. The treated cohort of CeD participants were consuming a GFD for at least 12 months. The active CeD cohort were recruited at diagnosis before starting a GFD and had positive transglutaminase-IgA or deamidated gliadin peptide-IgG and duodenal histology that subsequently confirmed CeD. Inclusion criteria for NCGS participants were self-report of gluten sensitivity, adherence to a GFD, documentation of prior exclusion of CeD based on negative CeD serology or small intestinal histology, or both, while eating gluten, or the presence of an HLA genotype not consistent with CeD. Inclusion criteria for nonceliac, non-gluten-sensitive controls were the absence of self-reported gluten sensitivity, regular gluten consumption, and documentation of prior exclusion of CeD confirmed by negative CeD serology or duodenal histology, or both, while eating gluten. Participants with autoimmune disease, such as type 1 diabetes mellitus and autoimmune thyroid disease, were diagnosed based on accepted clinical criteria. Exclusion criteria were the use of systemic immunosuppressant medication, pregnancy, or the presence of refractory CeD.

Procedures

Eligible participants attended a single visit where blood was collected for CeD serology (transglutaminase-IgA and deamidated gliadin peptide-IgG, Melbourne Health) and HLA-DQ2/DQ8 genotyping via single-nucleotide polymorphism tagging²¹ or polymerase chain reaction (Melbourne Pathology). Medical history and medications were recorded. A subset of participants with treated CeD, NCGS, and healthy controls on a GFD for at least 4 weeks underwent a single-dose open-label gluten challenge (10 g vital wheat gluten) for additional immune studies or to follow-up WBAIL-2 results discordant with their diagnosis. Serum for IL-2 assessment was collected at baseline and 4 hours after a gluten challenge (GCIL-2), while whole blood was collected at baseline (8 mL) and day 6 (8 mL) for WBAIL-2. In some cases, additional blood (50–300 mL) was collected at baseline and day 6 for tetramer analysis. A subset of treated CeD participants undertook a 3-day oral gluten

challenge (10 g/d) to enable cytokine capture assays on peripheral blood mononuclear cells (PBMCs) collected on day 6. Participants undergoing gluten challenges completed a diary (modified CeD patient-reported outcome)²² to record symptoms before, during, and after the challenge.

Cytokine Assays

The WBAIL-2 assay was adapted from a previously reported 96-well plate approach²⁰ to an individual tube format. The WBAIL-2 assay used in this study is a prototype version of the Novoleukin-C.G8 whole-blood assay, which is a proprietary commercial product developed and distributed by Novoviah Pharmaceuticals (Brisbane, Queensland, Australia). The stimulation tube contains a mixture of 8 synthetic immunoreactive gluten T-cell stimulatory peptides, including the well-characterized, overlapping DQ2.5-glia- α 1b, DQ2.5-glia- α 2, DQ2.5-glia- ω 1, and DQ2.5-glia- ω 2 epitopes included in the tetramers used in this study, as well as selected additional HLA-DQ2.5, HLA-DQ2.2, and HLA-DQ8-restricted epitopes derived from wheat α -, γ -, and ω -gliadins, and low-molecular-weight glutenin.⁸ The scrambled peptide negative control condition matched the amount and amino acid composition of the peptide mix in the stimulation tube. Synthetic peptides were >95% pure by high-performance liquid chromatography, and identity was confirmed by liquid chromatography-mass spectrometry (Mimotopes, Clayton, Victoria, Australia).

Whole-blood was collected in 4 mL heparinized tubes, and within 30 minutes, the peptide cocktail or a negative control (scrambled peptide) was added. Samples were incubated at 37°C for 24 hours $\pm$ 20 minutes. Plasma was then collected and stored at -80°C. IL-2 levels were measured using an electrochemiluminescence IL-2 S-PLEX kit (Meso Scale Diagnostics [MSD], Rockville, MD) following the manufacturer's instructions. A subset of samples was also assessed for interferon- γ (IFN- γ), IL17A, IL6, tumor necrosis factor- α , and IL10 using Proinflammatory and Cytokine V-PLEX assays. Plates were run on an SQ 120MM instrument (MSD), and mean cytokine levels from duplicate wells were analyzed using MSD Discovery Workbench software. The fold-change (peptide/control) and IL-2 concentration (peptide minus control) were calculated. The optimal diagnostic cutoff was determined post hoc using receiver operating characteristic curve analysis. Interassay and intra-assay variability were calculated (Supplementary Figure 1).

Tetramer Generation

The extracellular domain of HLA-DQ2.5 (HLA-DQA1*5:01 and HLA-DQB1*02:01) with the DQ2.5-hor-3a epitope or the DQ2.5-glia- α 1a, DQ2.5-glia- α 2, DQ2.5-glia- ω 1, DQ2.5-glia- ω 2 epitopes, or a class II-associated invariant chain peptide (CLIP; ATPLLMQALPMGA) control covalently linked to the N-terminus of the HLA-DQ2.5 β -chain via a flexible linker (GSGGSIEGRGSG) was produced and purified essentially as described.²³ Purified HLA-DQ2.5-peptide monomers were biotinylated on the C-terminal *Escherichia coli* BirA ligase recognition sequence (GLNDI-FEAQKIEWHE) of the HLA-DQ2.5 β -chain using BirA ligase buffer exchanged on a HiTrap Desalting column (Cytiva) to remove excess biotin, and complexed with streptavidin-phycoerythrin (PE) (BD Pharmingen) to form tetramers.

Tetramer Staining and Flow Cytometry

PBMCs were isolated from whole-blood samples by density gradient centrifugation (Leucosep, Interpath Services) and cryopreserved. PBMCs were incubated with 50 nmol/L dasatinib at 37°C for 30 minutes (Cell Signaling Technology), followed by an incubation with Fc receptor blocker (Miltenyi Biotec or Stem Cell Technologies). PBMCs were stained with HLA-DQ2.5-CLIP tetramer or HLA-DQ2.5:gluten tetramers conjugated to PE, representing the epitopes DQ2.5-glia- α 1a, DQ2.5-glia- α 2, DQ2.5-glia- ω 1, DQ2.5-glia- ω 2, and DQ2.5-hor3 (10 μ g/mL of each tetramer). Tetramer-stained PBMCs underwent tetramer-enrichment using the EasySep Release Human PE positive selection kit (Stem Cell Technologies) following the manufacturer's instructions (not the CLIP control). The CLIP negative control was used to establish a pre-enrichment number of CD4⁺ T cells. Cells were stained with Zombie UV live/dead marker (BioLegend) for 30 minutes at room temperature, followed by an antibody mix comprising anti-human CD3-Brilliant UltraViolet (BUV) 395 (clone: UCHT1; BD), CD4-Brilliant Violet (BV) 480 (clone: RPA-T4; BD), CD8-BV750 (clone: SK1; BioLegend), CD62L-fluorescein isothiocyanate (clone: DREG-56; BioLegend), CD14-allophycocyanin (APC)-H7 (clone: MDP9; BD Pharmingen), CD45RA-peridinin-chlorophyll protein (PerCP)-Cy5.5 (HI100; BD), Integrin β 7-BV421 (clone: FIB504; BD), and CD38-APC/Fire810 (clone: HB-7; BioLegend), prepared in Brilliant Stain Buffer (BD), for 30 minutes at 4°C. Cells were analyzed on a Cytex Aurora cytometer, and flow cytometry data were analyzed by FlowJo 10 software (BD Biosciences). CD4⁺ T cells were analyzed by flow cytometry for CD3⁺CD14⁻CD4⁺ HLA-DQ2.5:gluten tetramer⁺CD62L⁻CD45RA⁻Integrin β 7⁺, hereafter termed tetramer⁺ β 7⁺ T-effector memory (T_{EM}) cells. Expression of CD38⁺ tetramer⁺ β 7⁺ T_{EM} cells was assessed to determine activation status. The gating strategy is shown in Supplementary Figure 2. The number of tetramer⁺ β 7⁺ T_{EM} cells was normalized to 10⁶ CD4⁺ cells estimated from a pre-enriched sample.

Interleukin-2 Interferon- γ Capture Assay

Cytokine capture assay was performed on fresh or frozen PBMCs acquired from treated CeD patients 6 days after a 3-day oral gluten challenge (10 g vital wheat gluten daily). Briefly, as large numbers of PBMC were required, CD4⁺ T cells were isolated using EasySep Human CD4⁺ T Cell Isolation Kit (Stem Cell Technologies) then stained with HLA-DQ2.5-gluten tetramers and magnetically enriched, as above. Tetramer-enriched samples were incubated with 900,000 autologous PBMCs and glia- α 1/ α 2-, glia- ω 1/ ω 2- and hor-3a-containing 14 to 16 mer peptides (GL Biochem) at 15 μ g/mL equimolar in T-cell media (Iscove's Modified Dulbecco's Medium supplemented with 5% heat-inactivated pooled human serum [Australian Lifeblood], 1x GlutaMAX [Gibco], 1x nonessential amino acids [Gibco], and 50 μ mol/L 2-mercaptoethanol [Sigma-Aldrich]) for 4 hours at 37°C in 5% CO₂. A CLIP condition and unstimulated condition were included. Cytokine-secreting cells were identified using the IL-2-APC and IFN- γ -fluorescein isothiocyanate secretion assay kits (Miltenyi Biotec), following the manufacturer's instructions, with a shortened secretion period of 30 minutes. Cells were stained with Zombie UV Live/Dead marker (BioLegend) for 30 minutes at 4°C, followed by an antibody mix

comprising anti-human CD3-BUV395 (clone: UCHT1; BD), CD4-BV480 (clone: RPA-T4; BD), CD8-BV750 (clone: SK1; BioLegend), CD62L-BV510 (clone: DREG-56; BD), CD14-BUV661 (clone: M5E2), CD11c-BUV661 (clone: B-ly6; BD), CD19-BUV661 (clone: HIB19), CD56-BUV661 (clone: MY31; BD), CD45RA-PerCP-Cy5.5 (HI100, BD), integrin- β 7-BV421 (clone: FIB504; BD), CD38-APC/Fire810 (clone: HB-7; BioLegend), and CD60-BUV737 (clone: FN50; BD), prepared in Brilliant Stain Buffer (BD), for 30 minutes at 4°C. Cells were analyzed on a Cytex Aurora cytometer, and flow cytometry data were analyzed by FlowJo 10 software.

Interassay and Intra-assay Variability

The interassay and intra-assay variability of the WBAIL-2 and plasma IL-2 assessment (MSD S-PLEX) were evaluated (Supplementary Figure 1A). The interassay variability of the IL-2 whole-blood assay (WBA) was determined by testing samples on 2 occasions 1 to 2 weeks apart. The mean IL-2 fold change and pg/mL from triplicate WBA tests was used to compare across the 2 time points. The intra-assay variability was determined by triplicate WBAIL-2 tests performed from the same blood draw. Additionally, a different operator performed an additional triplicate set of WBAIL-2 tests to determine interoperator variability. The coefficient of variation (CV) was calculated as the ratio of the standard deviation to the mean, expressed as a percentage (%CV = standard deviation/mean $\times$ 100).

Statistical Analyses

Data were analyzed using GraphPad Prism 5 software (GraphPad Prism Software). Nonparametric tests were applied: the Mann-Whitney *U* test for unpaired data and Wilcoxon's signed rank test for paired data, with corrections for multiple comparisons.

Results

Participant Demographics

Between August 23, 2022, and December 18, 2024, 220 adults were screened for inclusion. Of these, 39 were excluded, and 181 were enrolled, comprising individuals with treated CeD ($n = 75$), active CeD ($n = 13$), NCGS ($n = 32$), and healthy controls ($n = 61$) (Table 1 and Supplementary Table 1). Women comprised 72% to 85% of the cohorts. Among the CeD participants, 85% (treated) and 92% (active) carried the celiac-susceptibility genotype HLA-DQ2.5 compared with 39% of healthy controls and 31% of NCGS participants.

The Whole-Blood Assay Measuring Interleukin-2 Release Is Sensitive and Specific for Celiac Disease

The diagnostic utility of WBAIL-2 was evaluated in a cohort of 88 adults with CeD (75 treated, 13 active) and 93 controls (32 NCGS, 61 healthy). IL-2 fold change and concentration were significantly higher in treated CeD (median, 7.77; 3.195 pg/mL) and active CeD (median, 3.73; 1.162 pg/mL) compared with healthy (median, 0.79; 0 pg/mL) and NCGS controls (median, 1.08; 0.011 pg/mL) (Figure 1A and B).

All CeD individuals were compared with all non-CeD individuals and an optimal cutoff of 1.99-fold change yielded 88% sensitivity (77 of 88) and 89% specificity (83 of 93; area under the curve [AUC], 0.92; 95% confidence interval [CI], 0.87–0.95; $P < .0001$) (Figure 1C; green). A secondary cutoff of 0.271 pg/mL achieved 90% sensitivity (79 of 88) and 83% specificity (77 of 93; AUC, 0.91; 95% CI, 0.86–0.96; $P < .0001$) (Figure 1D; green). Combining both cutoffs improved specificity to 95% (88 of 93) with 86% sensitivity (76 of 88) (Figures 1E and F). Treated CeD was also compared with NCGS to determine specificity in these 2 cohorts consuming a GFD. The 1.99-fold change cutoff yielded 81% specificity (26 of 32; AUC, 0.89; 95% CI, 0.82–0.95; $P < .0001$) (Figure 1C; blue), and the 0.271 pg/mL cutoff yielded 84% specificity (27 of 32; AUC, 0.90; 95% CI, 0.83–0.96; $P < .0001$) (Figure 1D; blue). The combined cutoff yielded 91% specificity (29 of 32) in participants consuming a GFD.

No significant differences in WBAIL-2 were observed across HLA types or zygosity (Figure 1G). However, sensitivity was lower (56%) in participants with the less common HLA-DQ8 genotype. Excluding this genotype from the analyses increased WBAIL-2 sensitivity to 90% (71 of 79). Using dual cutoff criteria in HLA-DQ2.5⁺ CeD participants, WBAIL-2 achieved 90% sensitivity (71 of 79) and 95% specificity (88 of 93).

Of 5 non-CeD participants with positive WBAIL-2 responses, 4 were first-degree relatives of a person with CeD (Table 2). First-degree relatives had a higher IL-2 fold change (median, 1.35) than nonrelatives (median, 0.74; $P = .0146$) (Figure 1H). Of these 5, 2 were healthy, and 3 reported NCGS. Three participants (1 healthy, 2 NCGS) agreed to a repeat WBAIL-2 assessment and underwent a single-dose gluten challenge after a 4-week GFD. The healthy participant and 1 NCGS participant again showed positive WBAIL-2 but had a negative GCIL-2. The other NCGS participant had a negative WBAIL-2 on retesting, suggesting an initial false positive.

Negative WBAIL-2 responses occurred in 12 CeD participants ($n = 10$ treated, $n = 2$ active) (Table 2). Among them were 4 HLA-DQ8⁺ individuals: 1 had no follow-up tests, 1 had prior positive IFN- γ ELISpot T-cell responses to rye and barley but not wheat, and 2 had high IL-2 concentrations in the WBA control tube. This pattern was also observed in 2 HLA-DQ2.5⁺ treated CeD nonresponders, who initially had substantially higher IL-2 concentrations in the control tube of their WBAIL-2 assays, but upon retesting, they had lower control tube IL-2 levels and their WBAIL-2 results were positive.

WBAIL-2 responses were initially negative in 2 treated CeD nonresponders, but these became positive on day 6 after a gluten challenge. One treated CeD participant had consistently negative WBAIL-2 results, including repeat WBAIL-2, postchallenge WBAIL-2, and GCIL-2. Both active CeD participants, after starting a GFD, showed positive WBAIL-2 responses when the test was repeated.

To evaluate the disease specificity of WBAIL-2, the CeD cohort was stratified by the presence ($n = 14$) or absence ($n = 61$) of concurrent autoimmunity. IL-2 fold change did not differ significantly between CeD participants with

Table 1. Demographics and Clinical Characteristics of Participants

Variable	Treated CeD (n = 75)	Active CeD (n = 13)	Healthy (n = 61)	NCGS (n = 32)
Female participants	64 (85)	10 (77)	44 (72)	27 (84)
Age, y	55 (20–76)	28 (18–52)	44 (19–71)	47 (24–75)
HLA-DQ				
2.5/2.5	7 (9)	2 (15)	3 (5)	0 (0)
2.5/x	35 (37)	7 (54)	14 (23)	9 (29)
2.5/2.2	14 (19)	2 (15)	2 (3)	0 (0)
2.5/7	3 (4)	0 (0)	1 (2)	0 (0)
2.5/8	7 (9)	1 (8)	4 (7)	1 (3)
2.2/2.2	0 (0)	0 (0)	7 (12)	0 (0)
2.2/x	1 (1)	0 (0)	1 (2)	6 (19)
2.2/8	0 (0)	0 (0)	0 (0)	1 (3)
8/8	2 (3)	0 (0)	0 (0)	0 (0)
8/7	1 (1)	1 (8)	0 (0)	1 (3)
8/x	5 (7)	0 (0)	9 (15)	3 (10)
7/x	0 (0)	0 (0)	6 (10)	7 (23)
x/x	0 (0)	0 (0)	13 (22)	3 (10)
Positive TTG and DGP serology	4 (5)	10 (77)	1 (2)	0 (0)
Positive TTG serology only	8 (11)	3 (23)	2 (3)	0 (0)
Positive DGP serology only	2 (3)	0 (0)	2 (3)	1 (3)
Years on a GFD, y	11 (1–52)	0	NA	7 (0.5–32)
Autoimmune disease other than CeD	14 (19)	0 (0)	10 (16)	6 (19)
Hashimoto's thyroiditis	5 (7)	0 (0)	3 (5)	4 (13)
Grave's disease	2 (3)	0 (0)	0 (0)	0 (0)
Rheumatoid arthritis	3 (4)	0 (0)	0 (0)	1 (3)
Type 1 diabetes mellitus	2 (3)	0 (0)	7 (11)	1 (3)
Dermatitis herpetiformis	1 (1)	0 (0)	0 (0)	0 (0)
Scleroderma	1 (1)	0 (0)	0 (0)	0 (0)
Addison's disease	1 (1)	0 (0)	0 (0)	0 (0)
Autoimmune hepatitis	1 (1)	0 (0)	0 (0)	0 (0)
Sjogren's syndrome	0 (0)	0 (0)	0 (0)	1 (3)
Lupus	0 (0)	0 (0)	0 (0)	1 (3)
GCIL-2	35	0	5	16
Tetramer analysis	10	0	0	4
Baseline WBAIL-2	75	13	61	32
Day 1 pregluten challenge WBAIL-2	31	0	1	12
Day 6 postgluten challenge WBAIL-2	31	0	1	11


NOTE. Data are presented as n (%), or median (range).

DGP, deamidated gliadin peptide antibody; NA, not applicable; NCGS, nonceliac gluten sensitivity; TTG, tissue transglutaminase antibody.

autoimmunity (median, 12.52) and those without (median, 7.66) (Figure 1I). Similarly, in the nonceliac cohort (n = 93), there was no significant difference between participants with (n = 16; median, 0.68) and without autoimmune disease (n = 77; median, 0.99) (Figure 1J). WBAIL-2 was positively correlated with years on a GFD ($r = 0.35$; 95% CI, 0.13–0.54; $P = .002$) (Figure 1K) and participant age ($r = 0.49$; 95% CI, 0.01–0.45; $P = .034$) (Figure 1L).

To determine whether the diagnostic utility of the WBAIL-2 could be improved, additional cytokines reported

to be increased by acute gluten challenge (IFN- γ , IL17A, IL6, tumor necrosis factor- α , and IL10)^{14,24} were measured in a subset of CeD (n = 43) and nonceliac (n = 26) participants. Samples were selected based on WBAIL-2 results either concordant (n = 58) or discordant (n = 11) with the diagnosis. Across cohorts, IL-2 correlated most strongly with IFN- γ ($r = 0.632$, $P = 2.2 \times 10^{-7}$) and IL-17A ($r = 0.595$, $P = 1.6 \times 10^{-7}$) (Supplementary Table 2). After IL-2, IFN- γ and IL17A showed the highest diagnostic performance (AUC, 0.69 and 0.64; $P = .018$ and $P = .091$,

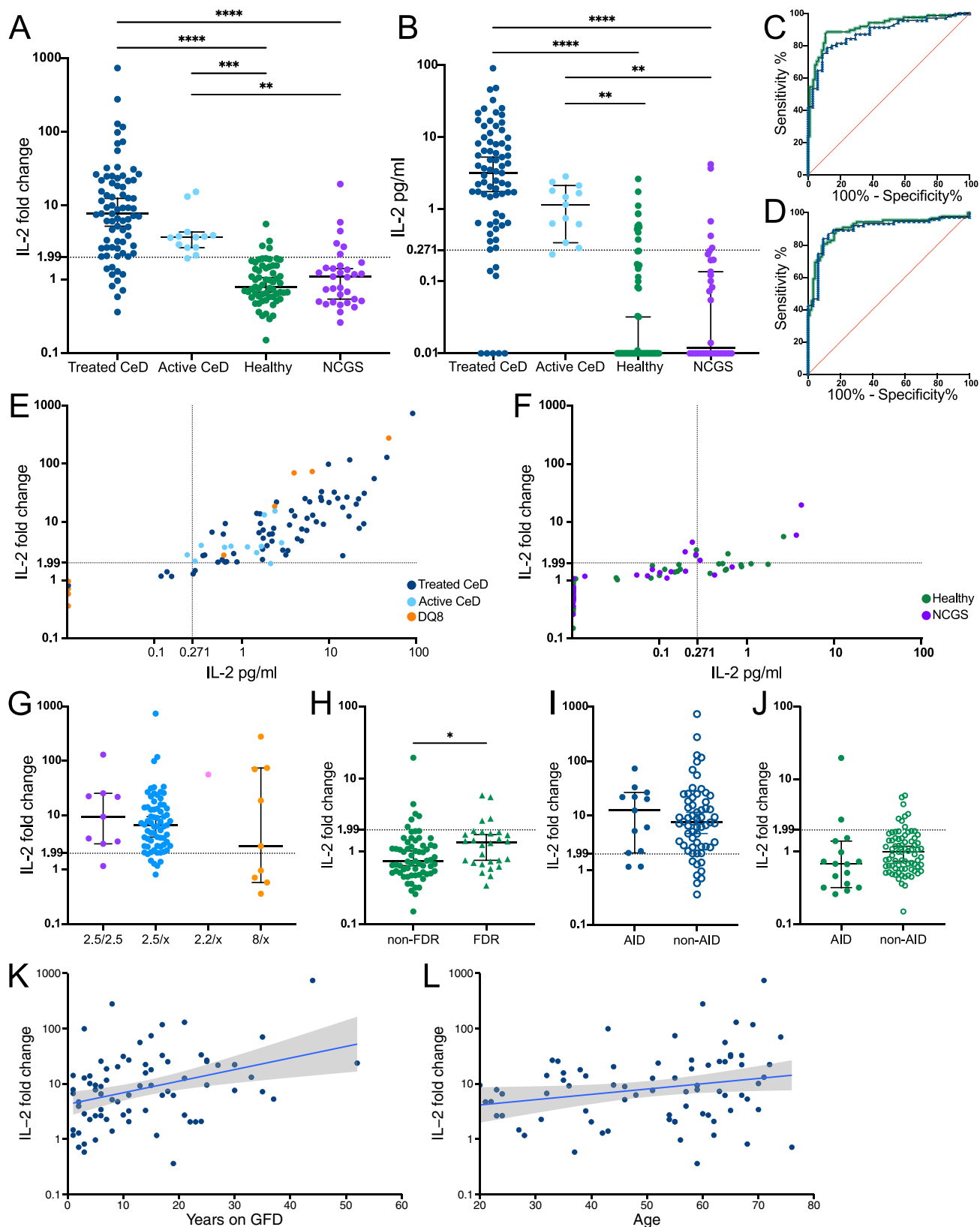

respectively). In participants with WBAIL-2 results consistent with their diagnosis, IFN- γ and IL17A had lower diagnostic accuracy than IL-2 (AUC, 0.84 and 0.75; $P = .0002$ and $P = .005$, respectively), and neither cytokine was positive in CeD participants with negative WBAIL-2 results. These findings confirm IL-2 as the optimal cytokine for CeD diagnostic assays. High IFN- γ responses in nonceliac participants with positive WBAIL-2 suggest a true biological response rather than technical issues or nonspecific IL-2 production.

The variability of WBAIL-2 and the IL-2 S-PLEX was assessed in 3 treated CeD participants and 3 nonceliac controls (Supplementary Figure 1A). For WBAIL-2, the mean CV was 24.1% to 29.7% for interassay variability, 28.9% to 46.3% for intra-assay variability, and 17.8% to 19.4% for interoperator variability (Supplementary Figure 1B-E). In the IL-2 S-PLEX, replicate wells had a mean CV of 2.5%, and interassay variability had a mean CV of 7.6% (Supplementary Figure 1F).

The Whole-Blood Assay Measuring Interleukin-2 Release Response Correlates With Serum Interleukin-2 After Oral Gluten Challenge

Serum IL-2 production after the oral gluten challenge (GCIL-2) was assessed to determine its diagnostic accuracy. GCIL-2 levels were significantly higher in CeD participants ($n = 35$; median, 16.1 pg/mL) compared with nonceliac participants ($n = 21$; median, 1.0 pg/mL; $P < .0001$) (Figure 2A). The optimal diagnostic cutoff for GCIL-2 was identified as a 2-fold increase from baseline to 4 hours postchallenge, achieving 97% sensitivity and 100% specificity for detecting CeD (AUC, 0.96; 95% CI, 0.9–1.0; $P < .0001$) (Figure 2B). Additionally, GCIL-2 levels were positively correlated with WBAIL-2 results ($r = 0.72$, 95% CI 0.56–0.83; $P < .0001$) (Figure 2C).

The Whole-Blood Assay Measuring Interleukin-2 Release Response Correlates With the Frequency of Activated Gluten-Specific CD4⁺ T Cells

The WBAIL-2 and GCIL-2 were compared with matched HLA-DQ2.5:gluten tetramer analysis. Higher tetramer⁺ $\beta 7^+$ T_{EM} cell frequencies in CeD participants ($n = 10$; median, 2.15) than in nonceliac participants ($n = 4$; median, 0.19;

$P = .024$) (Figure 2D) were shown. A cutoff of 0.59 cells per million CD4⁺ T cells achieved 100% sensitivity and 75% specificity for CeD (AUC, 0.90; 95% CI, 0.7–1.0; $P = .024$) (Figure 2E). Tetramer⁺ $\beta 7^+$ T_{EM} cell frequency correlated with GCIL-2 ($r = 0.74$, $P = .0035$) (Figure 2F) and with WBAIL-2 both with ($r = 0.79$; 95% CI, 0.58–0.90; $P < .0001$) (Figure 2G) and without ($r = 0.73$; 95% CI, 0.48–0.87; $P < .0001$) the gluten challenge (Figure 2H). A positive WBAIL-2 (≥ 2 -fold change over control) was detected with as few as 1 tetramer⁺ $\beta 7^+$ T_{EM} cell per million CD4⁺ T cells, equivalent to 1 cell in 4 mL of blood. WBAIL-2 also correlated with CD38⁺ tetramer⁺ $\beta 7^+$ T_{EM} cells ($r = 0.62$, $P = .0005$) (Figure 2I), representing activated gluten-specific T cells.

In Vivo Gluten Challenge Increases the Frequency of Tetramer⁺ $\beta 7^+$ T-Effector Memory Cells and the Whole-Blood Assay Measuring Interleukin-2 Release Response

WBAIL-2 and tetramer frequency were assessed at baseline (day 1) and day 6 after an oral gluten challenge in treated CeD. The gluten challenge expands gut-homing gluten-specific T cells in the blood by day 6.^{11,12,16,25,26} The WBAIL-2 increased by a median of 3-fold on day 6 compared with baseline, with IL-2 fold changes rising from a median of 7.8 on day 1 to 34.1 on day 6 ($P = .0013$) (Figure 3A). Two CeD participants with negative WBAIL-2 results on day 1 had positive responses on day 6, indicating gluten-induced expansion of gluten-specific T cells enhances in vitro IL-2 production. The oral gluten challenge increased absolute IL-2 concentrations in vitro by a median of 2.5-fold, from 1.5 pg/mL on day 1 to 7.4 pg/mL on day 6 ($P < .0001$) (Figure 3B). Similarly, the frequency of tetramer⁺ $\beta 7^+$ T_{EM} cells increased by a median of 2.4-fold, from 2.2 on day 1 to 8.7 on day 6 ($P = .0098$) (Figure 3C). CD38⁺ tetramer⁺ $\beta 7^+$ T_{EM} cells, indicative of recent gluten-induced activation,¹⁶ rose from 0.12 on day 1 to 6.8 on day 6 ($P = .004$) (Figure 3D). One participant showed a decrease in tetramer⁺ $\beta 7^+$ T_{EM} cells on day 6, consistent with their WBAIL-2 response. The percentage of tetramer⁺ $\beta 7^+$ T_{EM} cells expressing CD38 increased from 3.6% on day 1 to 76.6% on day 6 ($P = .0098$) (Figure 3E). These findings demonstrate that the WBAIL-2 reflects the expansion of circulating gluten-specific T cells after a gluten challenge.

Figure 1. WBAIL-2 diagnostic performance. (A) IL-2 fold change and (B) IL-2 concentration in treated CeD ($n = 75$), active CeD ($n = 13$), healthy controls ($n = 61$) and NCGS ($n = 32$). (C) Receiver operating characteristic curves are shown for IL-2 fold change and (D) concentration for treated CeD compared with NCGS and healthy (green) or compared with NCGS alone (blue). Performance of combined cutoffs using both IL-2 fold change and concentration in (E) CeD and (F) control cohorts. (G) IL-2 fold change in treated CeD based on HLA type and zygosity. (H) IL-2 fold change in non-CeD controls stratified by presence of a first-degree relative (FDR) with CeD. IL-2 fold change in (I) CeD and (J) non-CeD cohorts stratified by the presence or absence of autoimmune diseases (AID) other than CeD. IL-2 fold change correlation with (K) years consuming a GFD and (L) age. The cutoffs determined by the receiver operating characteristic curves are shown by the dashed lines. Negative IL-2 concentrations were set at 0.01 to allow graphing on a log₁₀ scale. Median and 95% CIs (shaded area) are shown. * $P < .05$, ** $P < .01$, *** $P < .001$, **** $P < .0001$.

Table 2. Participants With Discordant Whole-Blood Assay Measuring Interleukin-2 Release Responses

Participant	Baseline WBA	Repeat baseline WBA	Day 6 after challenge WBA	GCIL-2	HLA	Notes
Active CeD						
1	Negative	Positive ^a			2.5/8	Baseline WBA close to cutoff (1.93-fold change, 2.14 pg/mL)
2	Negative	Positive ^a			2.5/x	Baseline WBA close to cutoff (2.72-fold change, 0.24 pg/mL)
Treated CeD						
1	Negative				8/x	No follow-up tests due to death
2	Negative				8/x	Previously responded to rye and barley on ELISpot but not wheat
3	Negative	Positive	Positive	Positive	8/8	Baseline WBA high control IL-2 (0.75 pg/mL) compared with positive repeat WBA (0.22 pg/mL)
4	Negative	Negative			8/x	High control IL-2 in baseline (7.9 pg/mL) and repeat WBA (35.4 pg/mL)
5	Negative	Negative	Positive	Positive	2.5/x	
6	Negative	Negative	Positive	Positive	2.5/x	
7	Negative	Positive	Positive	Positive	2.5/x	Baseline WBA high control IL-2 (2.4 pg/mL) compared with positive repeat WBA (0.16 pg/mL)
8	Negative	Positive	Positive	Positive	2.5/2.2	Baseline WBA high control IL-2 (0.35 pg/mL) compared with positive repeat WBA (0.08 pg/mL)
9	Negative	Negative	Negative	Negative	2.5/2.5	In clinical trial on α -melanocyte-stimulating hormone for type 1 diabetes
10	Negative	Negative				No gluten challenge as pregnant
Healthy						
1	Positive	Positive	Positive	Negative	2.2/x	FDR, negative serology while eating gluten (2024)
2	Positive				2.5/x	FDR, negative serology and histology while eating gluten (2024)
NCGS						
1	Positive	Positive		Negative	7/x	FDR
2	Positive	Negative	Negative	Negative	8/x	FDR
3	Positive	Positive		Negative	7/x	

FDR, first degree relative; NCGS, nonceliac gluten sensitivity.

^aRepeat baseline WBAIL-2 done after 1 week on GFD for active CeD.

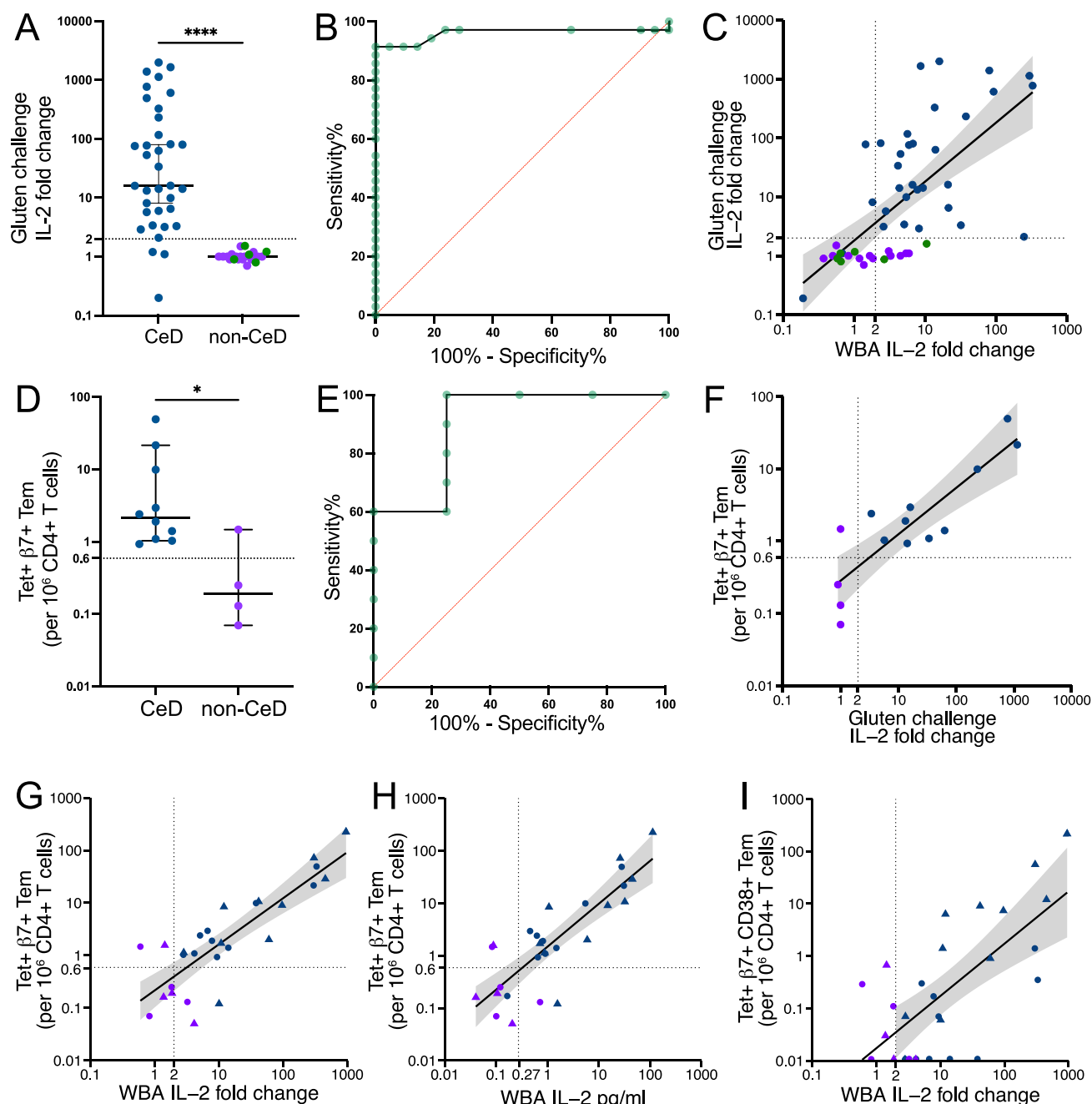


Figure 2. Comparison of WBAIL-2 with GCIL-2 and tetramer frequency. (A) GCIL-2 calculated as the fold change in IL-2 between baseline and 4 hours after the gluten challenge is shown for CeD ($n = 35$) and non-CeD cohorts ($n = 21$). (B) Receiver operating characteristic curve is shown for the GCIL-2. (C) GCIL-2 responses are plotted against WBAIL-2 responses (calculated as fold change between peptide and control tube). The cutoffs for a positive GCIL-2 and WBA (2-fold change) are indicated by the dotted lines. (D) The number of tetramer+ $\beta 7^+$ T_{EM} T cells per million CD4+ T cells are shown for CeD ($n = 10$) and non-CeD cohorts ($n = 4$). (E) Receiver operating characteristic curve is shown for the frequency of tetramer+ $\beta 7^+$ T_{EM} T cells per million CD4+ T cells. The frequency of tetramer+ $\beta 7^+$ T_{EM} T cells per million CD4+ T cells is plotted against (F) GCIL-2 fold change, (G) WBAIL-2 fold change, and (H) WBAIL-2 concentration (pg/mL). (I) Frequency of CD38+ tetramer+ $\beta 7^+$ T_{EM} T cells per million CD4+ T cells is plotted against the WBAIL-2 fold change. Treated CeD are shown in blue, NCGS in purple, and healthy participants in green, and baseline frequencies are indicated by circles, and frequencies on day 6 after the gluten challenge are indicated by triangles. Median and 95% CIs (shaded area) are shown. * $P < .05$, **** $P < .0001$.

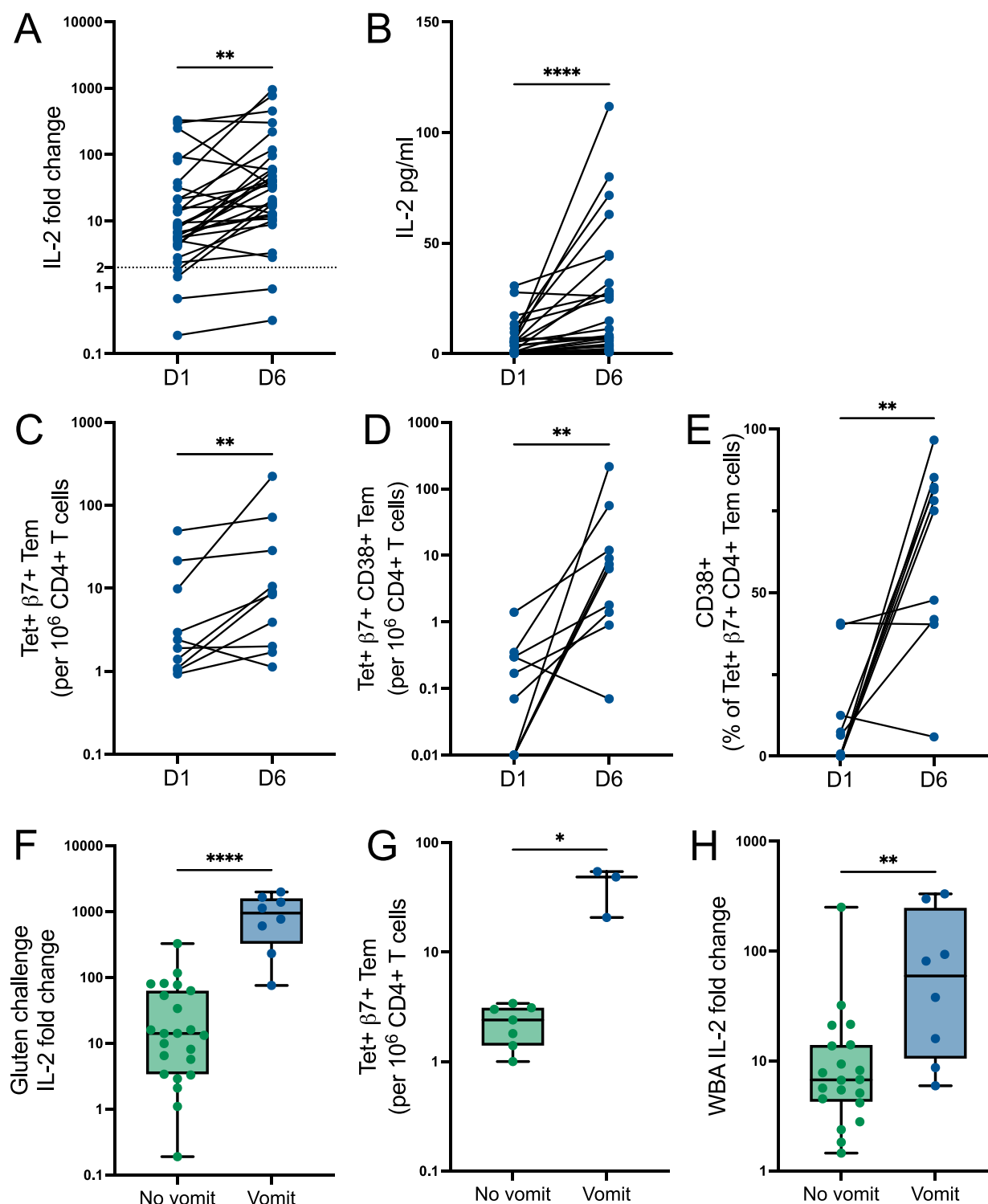


Figure 3. WBAIL-2 responses and tetramer frequencies reflect gluten challenge-induced T-cell expansion and symptoms in treated celiac disease. (A) WBAIL-2 fold change, (B) WBAIL-2 concentration (pg/mL), (C) frequency of tetramer⁺β7⁺T_{EM} T cells per million CD4⁺ T cells, (D) frequency of CD38⁺ tetramer⁺β7⁺T_{EM} T cells per million CD4⁺ T cells, and (E) percentage of tetramer⁺β7⁺T_{EM} T cells expressing CD38 are shown before (day 1; D1) and after (day 6; D6) the gluten challenge in treated CeD participants (n = 31 for WBAIL-2 and n = 10 for tetramer). Treated CeD participants were stratified based on whether or not they vomited after the oral gluten challenge. (F) GCIL-2 is shown, and both (G) frequency of tetramer⁺β7⁺T_{EM} T cells per million CD4⁺ T cells, and (H) WBA IL-2 performed before the gluten challenge are shown. Box and whisker plot: The boxes indicate the 25th percentile (bottom border), median (center line), and 75th percentile (top border), and the whiskers show the maximum and minimum ranges. **P* < .05, ***P* < .01, *****P* < .0001.

In Vivo and In Vitro Interleukin-2 Responses and Tetramer Frequency Correlate With Patient Symptoms

The relationship between gluten-induced symptoms and gluten-specific T-cell responses was assessed in 31 treated CeD participants. The most common symptoms after gluten exposure were tiredness (55%; mean severity, 4.5 of 10), nausea (55%; mean severity, 5.1 of 10), and headache (52%; mean severity, 3.7 of 10), with peak symptoms most commonly occurring at 3 to 4 hours (symptom heat map, [Supplementary Figure 3](#)). Severe symptoms (≥ 7 of 10) were reported by 9 participants (29%), 7 of whom experienced vomiting. Vomiting participants had significantly higher gluten-specific T-cell responses: GCIL-2 was 60-fold higher, with a median 949-fold change compared with 14-fold in nonvomiting participants ($P < .0001$) ([Figure 3F](#)), tetramer⁺ $\beta 7^+$ T_{EM} cell frequency was 20-fold higher, with a median of 48.5 compared with 2.4 ($P = .017$) ([Figure 3G](#)), and WBAIL-2 was 10-fold higher, with a median 60-fold change compared with 6.7-fold ($P = .004$) ([Figure 3H](#)). No consistent relationship was observed between other symptoms and T-cell responses, although notably, 1 nonvomiting participant with a >100 -fold WBAIL-2 increase reported severe tiredness ([Figure 3H](#)). These findings suggest that the magnitude of WBAIL-2 responses at baseline can identify treated CeD participants likely to experience severe symptoms after gluten exposure.

Gluten-Specific CD4⁺ T Cells Are Responsible for Gluten-Induced Interleukin-2 Production

To confirm acute IL-2 production is a biologically relevant marker of pathogenic gluten-specific T cells, a cytokine capture assay was performed using PBMCs from treated CeD participants after the gluten challenge ([Figure 4A](#)). From 125 to 205 million starting PBMCs, a median of 323 tetramer⁺ $\beta 7^+$ T_{EM} T cells was detected, with 64.7% (median, 209 cells) expressing CD69 ([Figure 4B](#)). The cytokine profile of CD69⁺tetramer⁺ $\beta 7^+$ T_{EM} T cells varied: 1 participant had more IL-2⁺IFN- γ ⁺ cells (38.9%) than IL-2⁺IFN- γ ⁺ cells (29.1%), and the other 2 had predominantly IL-2⁺IFN- γ ⁺ cells (61.7% and 28.7%) compared with IL-2⁺IFN- γ ⁺ cells (12.4% and 16.6%, respectively) ([Figure 4C](#)). All participants had few IL-2⁺IFN- γ ⁺ cells (median 4.97%). Higher frequencies of IL-2⁺ cells correlated with higher WBAIL-2 fold changes ([Figure 4D](#)). A similar trend was observed for IFN- γ , although 1 participant had a high frequency of IFN- γ ⁺ cells but low IFN- γ fold change, likely due to using fresh PBMCs instead of frozen, which may affect cytokine secretion. Shorter secretion periods favor IL-2 capture, whereas longer periods favor IFN- γ . No IL-2 was detected from naïve, central memory, or CD8⁺tetramer⁺ $\beta 7^+$ T cells (data not shown). These findings confirm that activated gluten-specific CD4⁺ T cells drive gluten-induced IL-2 production.

Discussion

Accurate blood-based T-cell diagnostics to support clinical practice are highly attractive, but none are used routinely outside of the QuantiFERON-Gold (Labcorp) assay for tuberculosis. CeD is an ideal candidate for T-cell diagnostics due to the central role of gluten-specific CD4⁺ T cells in its pathogenesis and the well-defined immunodominant gluten epitopes driving them.⁸ The field is increasingly shifting toward an immune-based diagnosis of CeD^{5,27,28} with increased acceptance of an immune (serology) non-biopsy specimen-based diagnosis when transglutaminase-IgA is 10 times the upper limit of normal in both children^{29,30} and, increasingly, adults.³¹ However, serology has limited positive-predictive value in average-risk patients and is uninformative for those avoiding gluten.³

Our findings show that the WBAIL-2 provides a direct measure of pathogenic gluten-specific T cells and offers high accuracy for diagnosing CeD, even in patients adhering to a strict GFD. The assay performed optimally in patients with the common HLA-DQ2.5 genotype, present in $>85\%$ of CeD cases.³² We confirm that IL-2 production in vitro correlates with serum IL-2 and the circulating frequency of gluten-specific T cells before and after gluten exposure. Strong WBAIL-2 responses predicted severe symptoms, such as vomiting, even before gluten ingestion. Importantly, we show that IL-2 is produced by tetramer⁺ gluten-specific T cells, highlighting the biological relevance of this diagnostic approach.

Detection of gluten-specific T cells is challenging due to their low frequency (~ 1 per 100,000 CD4⁺ T cells in treated CeD and negligible in individuals without CeD).^{10,16} To date, this has only been possible using HLA-DQ-gluten tetramer-based assays¹⁰ or prior expansion using a 3-day gluten challenge coupled with T-cell ELISpot or enzyme-linked immunosorbent assays.^{11,33} Although the tetramer-based approach incorporating 5 T-cell epitopes (HLA-DQ2.5-glia- $\alpha 1a$, DQ2.5-glia- $\alpha 2$, DQ2.5-glia- $\omega 1$, DQ2.5-glia- $\omega 2$, and DQ2.5-hor-3) has excellent sensitivity (97%) and specificity (95%) for CeD, the requirement for large blood volumes (>30 mL), knowledge of patient HLA genotype, technical expertise, and access to flow cytometry, makes this approach impractical for clinical translation.^{10,16} Further, approaches using a gluten challenge are less appealing to patients. The WBAIL-2 has similar sensitivity and specificity to these approaches but avoids the gluten challenge, is simple to perform, and requires limited blood (4 mL per tube).

In our community study, the WBAIL-2 demonstrated high sensitivity and specificity for the detection of CeD. The findings confirm the accuracy of the WBAIL-2 in a simple tube-based format. Variability (mean CV, 18%–46%) was similar to the QuantiFERON-Gold assay (mean CV, 13%–30%),³⁴ and likely reflects the challenges of detecting gluten-specific T cells that are circulating in low frequencies. Although GCIL-2 offered optimal diagnostic performance, the WBAIL-2 required only 8 mL of blood and no gluten challenge. We showed that priming with oral gluten

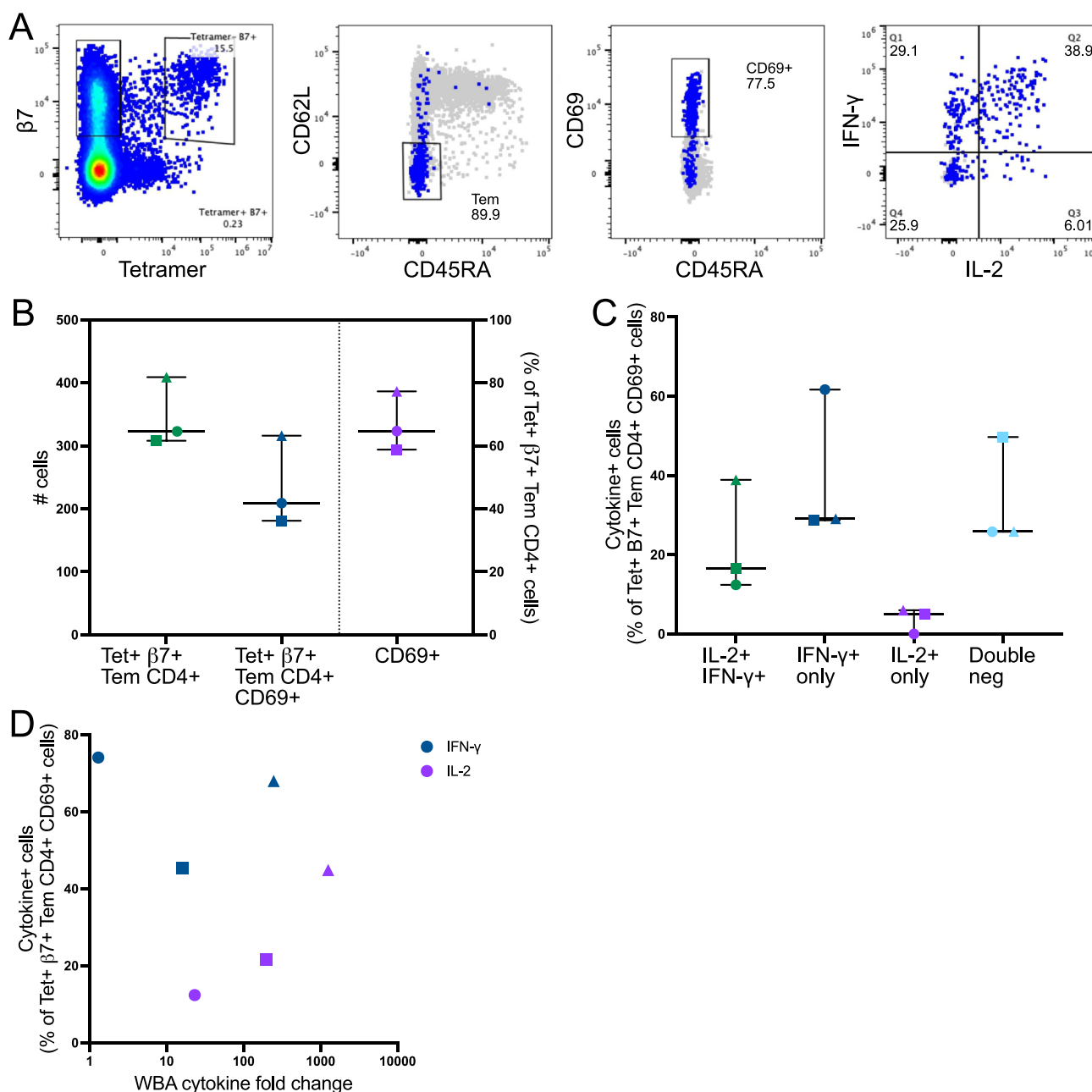


Figure 4. IL-2 and IFN- γ production by gluten-specific CD4⁺ T cells. (A) Gating strategy for IL-2- and IFN- γ -producing cells gated on tetramer⁺ $\beta 7$ ⁺ cells in blue and tetramer⁻ $\beta 7$ ⁺ in grey, CD45RA⁻CD62L⁻T_{EM}, and activated CD69⁺ cells, 1 representative of 3 similar experiments. (B) The number of tetramer⁺ $\beta 7$ ⁺T_{EM} cells and CD69⁺tetramer⁺ $\beta 7$ ⁺T_{EM} cells (left y-axis) and the frequency of tetramer⁺ $\beta 7$ ⁺T_{EM} T cells that express CD69 (right y-axis). (C) IL-2 and IFN- γ double positives, single positives, and double negatives as a frequency of CD69⁺ tetramer⁺ $\beta 7$ ⁺T_{EM} cells. (D) Frequency of CD69⁺ tetramer⁺ $\beta 7$ ⁺T_{EM} cells that express IFN- γ (blue) or IL-2 (purple) are plotted against the fold change in the respective cytokine from a WBA. Median and 95% CIs are shown. Each shape represents a different patient.

enhances the sensitivity of the WBAIL-2 by expanding the gluten-specific T-cell pool.

The lower sensitivity of the WBAIL-2 for HLA-DQ8⁺ CeD may have been impacted by the small sample size for this less common genotype, and another possibility is that HLA-DQ8-restricted responses in some CeD patients were below the assay's detection limit. Indeed, we previously showed that serum IL-2 responses after an oral gluten challenge in

HLA-DQ8⁺ CeD is lower than HLA-DQ2.5⁺ CeD, suggesting that HLA type can impact the T-cell response.¹⁸ Further, there may be other HLA-DQ8-restricted peptides that could enhance WBAIL-2 sensitivity, although detailed characterization of HLA-DQ8-restricted immunodominant epitopes is considerably more limited than in HLA-DQ2.5⁺ CeD.⁸

Prospective studies in larger populations are needed, both to better understand WBAIL-2 performance in HLA-

DQ8⁺ CeD and also inform the optimal testing algorithm for this assay. The WBAIL-2 could potentially be used as a confirmatory diagnostic in lieu of gastroscopy or as a first-line test in people consuming a GFD. The role of a single-dose gluten challenge to measure serum IL-2 or boost WBAIL-2 responses to clarify “equivocal” WBAIL-2 results also warrants investigation.

Confirming published work, we show a single-dose gluten challenge triggers a rise in serum IL-2 within hours, with high sensitivity (97%) and specificity (100%) for CeD, and an expansion of circulating CD38⁺ activated gluten-specific CD4⁺ T cells within days.^{14,16} Although WBAIL-2 closely correlates with GCIL-2 responses, some discordant results occur, with positive GCIL-2 but negative WBAIL-2, or vice versa. This may reflect differences in the assays: GCIL-2 measures systemic IL-2 responses to ingested wheat gluten, where a broad range of polyclonal gluten-specific T cells have been activated, whereas WBAIL-2 uses in vitro stimulation with a small number of deamidated gluten peptides, targeting a narrower T-cell repertoire in a low blood volume that may only contain a few gluten-specific T cells.

This study is the first to demonstrate that the gluten-specific T cell is responsible for gluten-induced IL-2 production, and it does this generally alongside IFN- γ production. These T cells fall into 3 functional subsets: IL-2⁺IFN- γ ⁻, IFN- γ ⁺IL-2⁻, and IFN- γ ⁺IL-2⁺, with the IFN- γ ⁺IL-2⁻ subset being most common in 2 of 3 participants. However, IL-2⁺IFN- γ ⁻ cells may be underestimated due to rapid IL-2 internalization,³⁵ and whether cytokine levels reflect a few highly active T cells or many less active ones remains unclear. Future work should examine how gluten-specific T cells shift between functional states over time and in relation to acute gluten exposure. These data confirm that gluten-induced IL-2 in the WBAIL-2 and GCIL-2 directly reflects the activation of pathogenic gluten-specific CD4⁺ T cells, confirming its biological relevance as a diagnostic biomarker in CeD.

We confirmed that higher GCIL-2 levels correlate with more severe symptoms^{14,36} and extended this to show that prechallenge WBAIL-2 responses and tetramer⁺ T-cell frequency predict gluten-induced vomiting. This suggests WBAIL-2 could predict symptom severity without an oral gluten challenge, aiding in selecting symptomatic patients for therapy trials where symptoms are key end points. Another advantage of the WBAIL-2 is monitoring the functional status of gluten-specific T cells without the need for a gluten challenge, which can be counterproductive and undesirable in the trial setting.

This study has several limitations. First, the single-center design and rigorous inclusion criteria limit the generalizability of the findings to diverse real-world patient populations, such as those taking systemic immunosuppression.

Second, the sample size for certain subgroups, particularly patients with less common non-HLA-DQ2.5 genotypes, was small, reducing the robustness of the conclusions drawn for these populations.

Third, although the WBAIL-2 demonstrated promising results, its reproducibility across different laboratories was

not assessed, which is critical for its adoption in routine clinical practice.

Fourth, the study also does not address the cost-effectiveness of implementing the WBAIL-2 compared with current diagnostic methods, an important consideration for real-world applicability. Studies assessing IL-2-based T-cell diagnostics in children are also needed, especially as we showed WBAIL-2 responses correlated positively with older participant age and previously showed serum IL-2 responses to oral gluten may be lower in younger patients with CeD.¹⁸ Future prospective validation studies with multiple centers and larger patient numbers, including those with non-HLA-DQ2.5 genotypes, will provide important data to validate the diagnostic performance of the WBAIL-2 in different clinical settings.

Using the WBAIL-2 approach, our findings highlight the potential for practical, blood-based T-cell diagnostics to extend beyond infectious diseases into autoimmunity. The ultrasensitive electrochemiluminescent platform is well-suited to clinical laboratories due to its compact design, minimal setup requirements, and automation that reduces the need for extensive technical training. Its high sensitivity, broad dynamic range, and capability to analyze various sample types, including serum and plasma, make it a practical choice for diagnostic purposes. Such a platform could be adapted for a variety of T cell-driven diseases and conditions, including type 1 diabetes, malignancies, transplant immunology, and other infectious diseases, whenever circulating antigen-specific T cells are likely to be present, even at low frequencies. Future studies should explore how whole-blood cytokine release assays, combined with ultrasensitive cytokine detection, can further advance T-cell diagnostics in these contexts, improving early detection, disease monitoring, and therapeutic stratification.

Conclusion

Our study indicates that gluten-stimulated IL-2 release assays are practical and correlate with the presence of a pathogenic gluten-specific CD4⁺ T-cell response in CeD. Availability of complementary in vivo and in vitro IL-2 release assays in CeD could address the diagnostic needs of clinicians managing patients with an uncertain diagnosis while avoiding a prolonged gluten rechallenge, children and adults unsuitable or unwilling to undergo endoscopy for histologic diagnosis, and provides a biomarker for stratifying disease severity by the “strength” of antigluten immunity. In immune conditions with T cells specific for well-defined antigens, in-tube WBAIL-2 release assays could be feasible and facilitate personalized diagnosis and therapy previously not possible in clinical care.

Supplementary Material

Note: To access the supplementary material accompanying this article, visit the online version of *Gastroenterology* at www.gastrojournal.org, and at <https://doi.org/10.1053/j.gastro.2025.05.022>.

References

- Iversen R, Sollid LM. The immunobiology and pathogenesis of celiac disease. *Annu Rev Pathol* 2023;18:47–70.
- Singh P, Arora A, Strand TA, et al. Global prevalence of celiac disease: systematic review and meta-analysis. *Clin Gastroenterol Hepatol* 2018;16:823–836.e2.
- Bibbins-Domingo K, Grossman DC, et al; US Preventive Services Task Force. Screening for celiac disease: US Preventive Services Task Force Recommendation Statement. *JAMA* 2017;317:1252–1257.
- Kvamme JM, Sorbye S, Florholmen J, et al. Population-based screening for celiac disease reveals that the majority of patients are undiagnosed and improve on a gluten-free diet. *Sci Rep* 2022;12:12647.
- Tye-Din JA. Evolution in coeliac disease diagnosis and management. *JGH Open* 2024;8:e13107.
- Lalvani A, Pathan AA, McShane H, et al. Rapid detection of *Mycobacterium tuberculosis* infection by enumeration of antigen-specific T cells. *Am J Respir Crit Care Med* 2001;163:824–828.
- Mori T, Sakatani M, Yamagishi F, et al. Specific detection of tuberculosis infection: an interferon-gamma-based assay using new antigens. *Am J Respir Crit Care Med* 2004;170:59–64.
- Sollid LM, Tye-Din JA, Qiao SW, et al. Update 2020: nomenclature and listing of celiac disease-relevant gluten epitopes recognized by CD4(+) T cells. *Immunogenetics* 2020;72:85–88.
- Christophersen A, Raki M, Bergseng E, et al. Tetramer-visualized gluten-specific CD4+ T cells in blood as a potential diagnostic marker for coeliac disease without oral gluten challenge. *United European Gastroenterol J* 2014;2:268–278.
- Sarna VK, Lundin KEA, Morkrid L, et al. HLA-DQ-gluten tetramer blood test accurately identifies patients with and without celiac disease in absence of gluten consumption. *Gastroenterology* 2018;154:886–896.e6.
- Anderson RP, Degano P, Godkin AJ, et al. In vivo antigen challenge in celiac disease identifies a single transglutaminase-modified peptide as the dominant A-gliadin T-cell epitope. *Nat Med* 2000;6:337–342.
- Tye-Din JA, Stewart JA, Dromey JA, et al. Comprehensive, quantitative mapping of T cell epitopes in gluten in celiac disease. *Sci Transl Med* 2010;2:1–14.
- Goel G, King T, Daveson AJ, et al. Epitope-specific immunotherapy targeting CD4-positive T cells in coeliac disease: two randomised, double-blind, placebo-controlled phase 1 studies. *Lancet Gastroenterol Hepatol* 2017;2:479–493.
- Goel G, Tye-Din JA, Qiao SW, et al. Cytokine release and gastrointestinal symptoms after gluten challenge in celiac disease. *Sci Adv* 2019;5:eaaw7756.
- Leonard MM, Silvester JA, Leffler D, et al. Evaluating responses to gluten challenge: a randomized, double-blind, 2-dose gluten challenge trial. *Gastroenterology* 2021;160:720–733.e8.
- Zuhlke S, Risnes LF, Dahal-Koirala S, et al. CD38 expression on gluten-specific T cells is a robust marker of gluten re-exposure in coeliac disease. *United European Gastroenterol J* 2019;7:1337–1344.
- Cartee AK, Choung RS, King KS, et al. Plasma IL-2 and symptoms response after acute gluten exposure in subjects with celiac disease or nonceliac gluten sensitivity. *Am J Gastroenterol* 2022;117:319–326.
- Tye-Din JA, Daveson AJM, Goldstein KE, et al. Patient factors influencing acute gluten reactions and cytokine release in treated coeliac disease. *BMC Med* 2020;18:362.
- Tye-Din JA, Skodje GI, Sarna VK, et al. Cytokine release after gluten ingestion differentiates coeliac disease from self-reported gluten sensitivity. *United European Gastroenterol J* 2020;8:108–118.
- Anderson RP, Goel G, Hardy MY, et al. Whole blood interleukin-2 release test to detect and characterize rare circulating gluten-specific T cell responses in coeliac disease. *Clin Exp Immunol* 2021;204:321–334.
- Hardy MY, Ontiveros N, Varney MD, et al. Resolving incomplete single nucleotide polymorphism tagging of HLA-DQ2.2 for coeliac disease genotyping using digital droplet PCR. *HLA* 2018;91:280–288.
- Leffler DA, Kelly CP, Green PH, et al. Larazotide acetate for persistent symptoms of celiac disease despite a gluten-free diet: a randomized controlled trial. *Gastroenterology* 2015;148:1311–1319.e6.
- Petersen J, Montserrat V, Mujico JR, et al. T-cell receptor recognition of HLA-DQ2-gliadin complexes associated with celiac disease. *Nat Struct Mol Biol* 2014;21:480–488.
- Goel G, Daveson AJM, Hooi CE, et al. Serum cytokines elevated during gluten-mediated cytokine release in coeliac disease. *Clin Exp Immunol* 2020;199:68–78.
- Raki M, Fallang LE, Brottveit M, et al. Tetramer visualization of gut-homing gluten-specific T cells in the peripheral blood of celiac disease patients. *Proc Natl Acad Sci U S A* 2007;104:2831–2836.
- Anderson RP, van Heel DA, Tye-Din JA, et al. T cells in peripheral blood after gluten challenge in coeliac disease. *Gut* 2005;54:1217–1223.
- Anderson RP. Emergence of an adaptive immune paradigm to explain celiac disease: a perspective on new evidence and implications for future interventions and diagnosis. *Expert Rev Clin Immunol* 2022;18:75–91.
- Lebwohl B, Ma C, Lagana SM, et al. Standardizing randomized controlled trials in celiac disease: an international multidisciplinary appropriateness study. *Gastroenterology* 2024;166:88–102.
- Husby S, Koletzko S, Korponay-Szabo I, et al. European Society Paediatric Gastroenterology, Hepatology and nutrition guidelines for diagnosing coeliac disease 2020. *J Pediatr Gastroenterol Nutr* 2020;70:141–156.
- Hoyle A, Gillett P, Gillett HR, et al. No-biopsy strategy for coeliac disease is applicable in adult patients: a 'real-world' Scottish experience. *Frontline Gastroenterol* 2023;14:97–102.
- Shiha MG, Nandi N, Raju SA, et al. Accuracy of the no-biopsy approach for the diagnosis of celiac disease in adults: a systematic review and meta-analysis. *Gastroenterology* 2024;166:620–630.
- Karell K, Louka AS, Moodie SJ, et al. HLA types in celiac disease patients not carrying the DQA1*05-DQB1*02

- (DQ2) heterodimer: results from the European genetics cluster on celiac disease. *Hum Immunol* 2003; 64:469–477.
33. Ontiveros N, Tye-Din JA, Hardy MY, et al. Ex-vivo whole blood secretion of interferon (IFN)-gamma and IFN-gamma-inducible protein-10 measured by enzyme-linked immunosorbent assay are as sensitive as IFN-gamma enzyme-linked immunospot for the detection of gluten-reactive T cells in human leucocyte antigen (HLA)-DQ2.5(+)-associated coeliac disease. *Clin Exp Immunol* 2014; 175:305–315.
 34. Tagmouti S, Slater M, Benedetti A, et al. Reproducibility of interferon gamma (IFN-gamma) release Assays. A systematic review. *Ann Am Thorac Soc* 2014;11:1267–1276.
 35. Deng N, Mosmann TR. Optimization of the cytokine secretion assay for human IL-2 in single and combination assays. *Cytometry A* 2015;87:777–783.
 36. Daveson AJM, Tye-Din JA, Goel G, et al. Masked bolus gluten challenge low in FODMAPs implicates nausea and vomiting as key symptoms associated with immune activation in treated coeliac disease. *Aliment Pharmacol Ther* 2020;51:244–252.

Received March 23, 2025. Accepted May 2, 2025.

Correspondence

Address correspondence to: Jason Tye-Din, MBBS, FRACP, PhD, Walter and Eliza Hall Institute, 1G Royal Parade, Parkville, Victoria 3052, Australia. e-mail: tyedin@wehi.edu.au.

Acknowledgments

The authors thank all of the volunteers who participated in the study. They thank Sinead Browne for help with study recruitment and sample collection, Coeliac Australia for funding support, and Novoviah Pharmaceuticals, who provided the proprietary tubes and peptides for this study.

CRedit Authorship Contributions

Olivia G. Moscatelli, BSc (Conceptualization: Equal; Formal analysis: Lead; Investigation: Lead; Visualization: Lead; Writing – original draft: Equal; Writing – review & editing: Supporting)

Amy K. Russell, BSc (Investigation: Supporting; Writing – review & editing: Supporting)

Lee M. Henneken, BN, BBiolSci (Investigation: Supporting; Project administration: Lead; Writing – review & editing: Supporting)

Linda Fothergill, PhD (Investigation: Supporting; Writing – review & editing: Supporting)

Allan Motyer, PhD (Formal analysis: Supporting; Visualization: Supporting; Writing – review & editing: Supporting)

Hugh Reid, PhD (Resources: Supporting; Writing – review & editing: Supporting)

Jamie Rossjohn, PhD (Resources: Lead; Writing – review & editing: Supporting)

Vanessa Bryant, PhD (Supervision: Equal; Writing – review & editing: Supporting)

Robert P. Anderson, MBChB, FRACP, PhD, AGAF (Conceptualization: Equal; Funding acquisition: Supporting; Writing – review & editing: Supporting)

Melinda Y. Hardy, PhD (Investigation: Supporting; Supervision: Equal; Writing – review & editing: Supporting)

Jason Allan Tye-Din, MBBS, FRACP, PhD (Conceptualization: Equal; Funding acquisition: Lead; Supervision: Equal; Writing – original draft: Equal; Writing – review & editing: Supporting)

Conflicts of interest


These authors disclose the following: Melinda Y. Hardy is a consultant for Takeda. Robert P. Anderson is an advisor to Immunic, Takeda Pharmaceutical, AMYRA Biotech, DBV Technologies, Barinthus Biotherapeutics, EVOQ Therapeutics, Pectigen ActoBio, Allero Therapeutics, Kanyos Bio, GSK, Treg Therapeutics, Bionix Therapeutics, AdAlta, Universal Cells, and Bristol-Myers Squibb Australia; is a director and shareholder of Novoviah Pharmaceuticals; and is inventor of patents relating to diagnosis and treatment of celiac disease. Jason Allan Tye-Din has privately or via his institution been a consultant or advisory board member for Anantara, Anokion, Barinthus Biotherapeutics, Chugai Pharmaceuticals, DBV Technologies, Dr. Falk, EVOQ Therapeutics, Equillium, Forte Biosciences, IM Therapeutics, Janssen, Kallyope, Mozart Therapeutics, Takeda, TEVA, and Topas, has received research funding from Barinthus Biotherapeutics, Chugai Pharmaceuticals, Codexis, DBV Technologies, Kallyope, Novoviah Pharmaceuticals, Topas, and Tillotts Pharmaceuticals, and is an inventor on patents relating to the use of gluten peptides in celiac disease diagnosis and treatment. The remaining authors disclose no conflicts.

Funding


Funding support for this study was provided by Coeliac Australia, Beck Family Foundation, Butterfield Family, and the Veith Foundation. Melinda Y. Hardy was supported by National Health and Medical Research Council Ideas Grant APP2012225 and Jason Allan Tye-Din and Jamie Rossjohn by National Health and Medical Research Council Investigator Grants APP1176553 and APP2008981. The content of this research is solely the responsibility of the authors; the funding bodies and Novoviah Pharmaceuticals played no role in the execution or interpretation of the research.

Data Availability


The data that support the findings of this study are available from the corresponding author upon reasonable request.

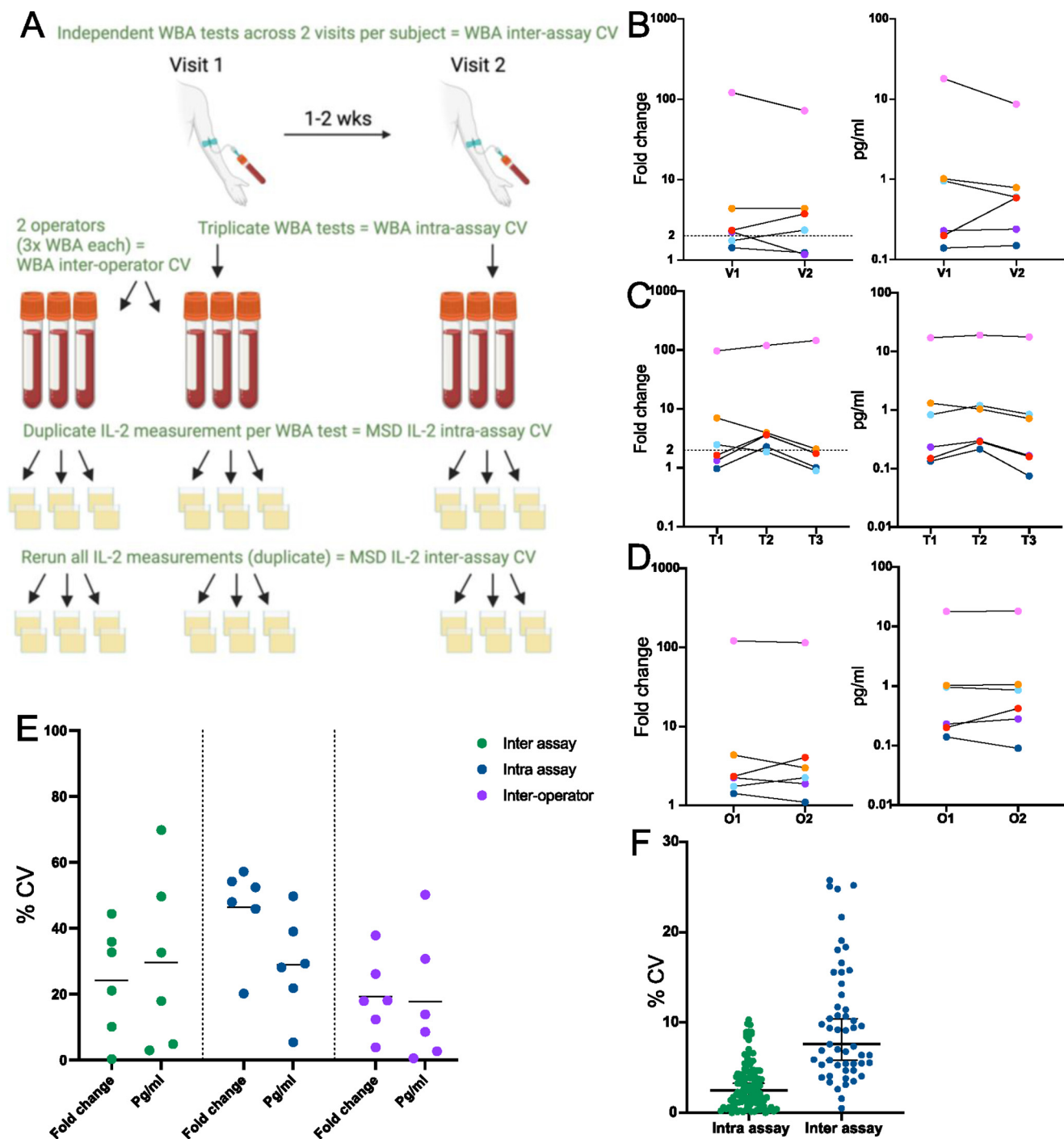

Supplementary Figure 1. Evaluation of interassay and intra-assay variability of the WBAIL-2 and the MSD IL-2 S-PLEX. (A) Experimental design for assessment of interassay and intra-assay and interoperator variability. IL-2 fold change (*left*) and concentration (pg/mL; *right*) are shown for (B) interassay replicates at visit 1 and 2, (C) intra-assay replicates for tests 1, 2 and 3, and (D) interoperator replicates for operators 1 and 2 for the WBAIL-2. (E) Percentage CVs are shown for fold change and pg/mL for interassay, intra-assay, and interoperator variability of the WBAIL-2. (F) Percentage CVs are shown for intra-assay and interassay variability of the MSD IL-2 S-PLEX.

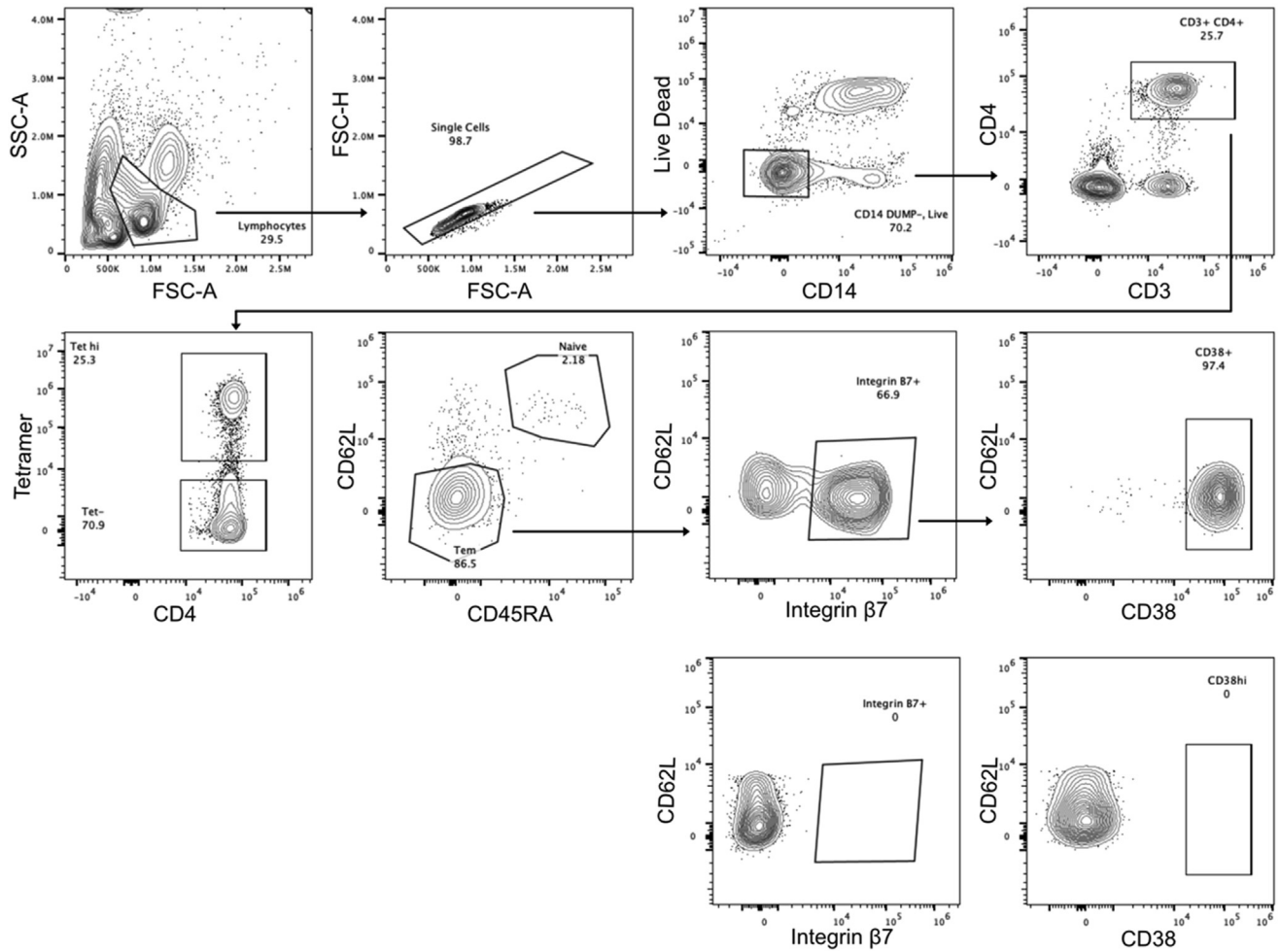

Supplementary Figure 2. Tetramer gating strategy. Example gating strategy of the tetramer-enriched sample. Shown is 1 representative treated CeD patient on day 6 after the gluten challenge. Firstly, cleanup gates were applied to gate lymphocytes, single cells, live CD14⁻ cells and CD3⁺CD4⁺ cells (*top row*). FSC-A, forward side scatter area; SSC-A, side scatter area. Next, tetramer⁺ cells were gated on, followed by effector memory cells identified as CD45RA⁻CD62L⁻, then integrin $\beta 7$ ⁺ cells, and lastly, CD38^{hi} cells (*middle row*). Integrin $\beta 7$ ⁺ and CD38⁺ gates are shown for respective fluorescence minus one samples (*bottom row*).

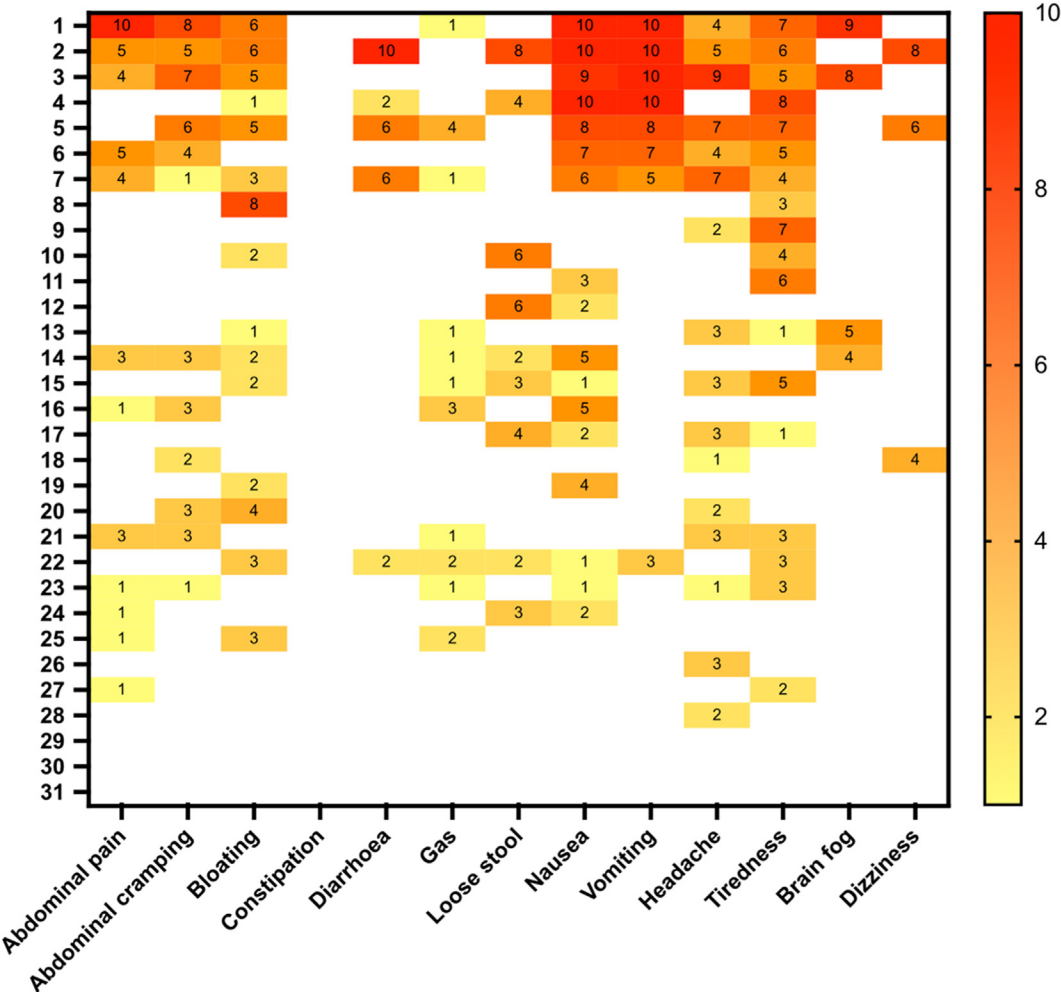

Supplementary Figure 3. Heat map of gluten-induced symptoms. Heat map of patient-reported symptom scores after gluten challenge in treated CeD participants (n = 31). Each row represents 1 participant. Maximum reported severity scores are shown, ranging from 0 (no symptoms; white) to 10 (severe symptoms; red).

Supplementary Table 1. Participant Screening, Inclusion and Exclusion Numbers

Variable	Treated CeD (n)	Active CeD (n)	Healthy (n)	NCGS (n)	Total (n)
Screened	96	13	66	45	220
Included	75	13	61	32	181
Excluded	21	0	5	13	39
Reason for exclusion					
Use of immunosuppressants	5	0	2	1	8
Uncertain diagnosis ^a	12	0	1	8	21
Laboratory technical issues	1	0	2	1	4
Nonadherent with GFD	3	0	0	3	6

NCGS, nonceliac gluten sensitivity.

^aUncertain diagnosis includes insufficient documentation to confirm the diagnosis, diagnosis on biopsy or serology alone and serology/biopsy specimen testing performed while the patient was consuming a GFD.

Supplementary Table 2. Additional Cytokine Analysis Performed on the Whole-Blood Assay Plasma

Variable	Correlation with IL-2				Diagnostic performance					
	pg/mL		Fold change		All participants (N = 69)			Participants with WBAIL-2 consistent with diagnosis (n = 58)		
	Spearman's <i>r</i>	<i>P</i> value	Spearman's <i>r</i>	<i>P</i> value	AUC	95% CI	<i>P</i> value	AUC	95% CI	<i>P</i> value
TNF- α	0.257	.033	0.348	.003	0.58	0.44-0.71	.287	0.64	0.50-0.79	.066
IL6	0.293	.015	0.261	.030	0.62	0.49-0.76	.090	0.66	0.52-0.80	.041
IL10	0.218	.072	0.225	.063	0.55	0.41-0.69	.492	0.61	0.46-0.76	.157
IFN- γ	0.684	9.66×10^{-11}	0.706	1.27×10^{-11}	0.75	0.64-0.86	.0005	0.87	0.77-0.96	<.0001
IL17A	0.683	1.02×10^{-10}	0.595	6.80×10^{-8}	0.73	0.61-0.85	.0016	0.82	0.71-0.93	<.0001